

Roche Applied Science
LightCycler[®] 480
LIMS Module
programming reference manual



Roche Applied Science
LightCycler[®]480
LIMS Module
programming reference manual

LIMS Software Version 1.5

*This is the paper version of the LIMS manual for the
Roche LightCycler[®]480*

*An electronic version in the HTML format is also
maintained and available from the installation CD
delivered with the instruments, or available on request to
a Roche representative.*

Table of Contents

Part I Introduction	8
1 Target audience	8
2 Document Revision History	8
3 Contacts, Trademarks	9
4 Glossary	10
Part II LIMSCientLib Reference	11
1 General information and debugging	12
General information	12
Install LIMS server on another PC	16
Multiple LIMS users	18
2 Restrictions and known issues	19
Environment	19
Other known issues	19
3 Interfaces	20
ILIMSConnection	21
LoggedIn Property.....	22
Query Property.....	23
ExperimentInfo Property.....	24
Instrument Property.....	25
Host Property	26
Port Property.....	27
Login Method.....	28
Logout Method.....	29
ILIMSExperimentInfo	30
GetStatus Method.....	31
GetCompletedExperimentSummary Method	33
ExportExperiment Method.....	35
ILIMSIInstrument	36
Reserve Method.....	37
Unreserve Method.....	38
Open Method.....	39
Close Method.....	41
StartExperiment method.....	42
GetStatus Method.....	44
AbortExperiment Method.....	46
GetContainerBarcode Method.....	47
OpenAndWait Method.....	48
CloseAndWait Method	49
GetContainerSensor.....	50
SetContainerSensor.....	52
ILIMSOperationResult	53
Successful Property.....	54

	ServerError Property	55
	Message Property	55
	ErrorNumber Property	56
	DateTime Property	56
	UserMessage Property	57
	ILIMSSampleDefinition	58
	SampleCount Property	59
	ExpDataCount Property	60
	AddSample Method	61
	Clear Method	63
	DeleteSample Method	64
	GetSample Method	65
	AddExpData Method	67
	GetExpData Method	68
	DeleteExpData Method	69
	ClearExpData Method	70
	ILIMSSampleInfo	71
	Position Property	72
	Name Property	73
	ID Property	74
	Comment Property	75
	ReplicatePosition Property *****	76
	ExpDataCount Property	77
	AddExpData Method	77
	GetExpData Method	78
	DeleteExpData Method	78
	ClearExpData Method	79
	ILIMSQuery	80
	ObjectType Property	81
	Name Property	82
	Owner Property	82
	FromDate property	83
	ToDate Property	84
	QueryDateType property	85
	ExecuteQuery Method	86
	ILIMSQueryResult	87
	Count property	88
	GetResultData method	89
	ILIMSQueryResultData	90
	Name Property	91
	ObjectType Property	91
	CreationDate Property	92
	ModificationDate Property	92
	Path Property	93
4	Enumerated Types	94
	LIMSQueryDateType	94
5	CoClasses	95
	LIMSConnection	96
	LIMSExperimentInfo	96
	LIMSInstrument	96
	LIMSOperationResult	96
	LIMSSampleDefinition	96

	LIMSSampleInfo	96
	LIMSQuery	96
	LIMSQueryResult	96
	LIMSQueryResultData	97
Part III	XML output	98
1	Additional details, requirements	98
2	Output file schema	100
3	XML-File Elements	101
	Experiment	101
	Run.....	101
	Protocol.....	102
	Programs.....	102
	HTCRunProgram nodes (Emlist)	102
	HTCRunSegment.....	103
	AnalysisModes.....	103
	DetectionFormat.....	103
	HTCDetectionFormat	103
	BlockType	103
	Acquisition.....	104
	Acquisition.....	104
	TempLog.....	104
	Analyses.....	105
	Analysis.....	105
	RelQuantGroup.....	105
	ReferenceGroup/TargetGroup	105
	RelQuantGroupData.....	105
	RelQuantSamples.....	105
	Pairings.....	105
	non RelQuant Analysis.....	107
	AnalysisSamples (Abs Quant/Fit Points).....	107
	AnalysisSample (Genotyping).....	107
	AnalysisSample (TM Calling).....	107
	AnalysisSample (Color Compensation).....	108
	Subsets.....	109
	Samples.....	109
Index		111

1 Introduction

1.1 Target audience

The intended audience for this manual is programmers possessing a good knowledge with Windows API programming.

The present manual will help such programmers to understand the basics of the LightCycler 480 LIMS interface and build efficient applications to control the instrument and collect the results of experiments.

Code examples are presented in this manual, to assist in the learning process. The functionality of those examples has been tested, however, as they have been edited to accommodate the manual constraints, there is no guaranty that they will directly work by simply using a "cut and paste" process into a practical application. In some instances, a slight adaptation might be required.

1.2 Document Revision History

Document Version	Revision Date
1.0	January 2008

Copyright 2007, Roche Diagnostics Ltd. All rights reserved.

Information in this document is subject to change without notice. No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Roche Diagnostics GmbH.

Questions or comments regarding the contents of this manual can be directed to the address below or to your Roche representative.

Roche Diagnostics Ltd.
Roche Applied Science
Global Service Support
Forrenstrasse
6343 Rotkreuz, Switzerland

Every effort has been made to ensure that all the information contained in this manual is correct at the time of printing.
However, Roche Diagnostics Ltd. reserves the right to make any changes necessary without notice as part of ongoing product development.

1.3 Contacts, Trademarks

Contact Addresses



Manufacturer

Roche Diagnostics Ltd.
Forrenstrasse
CH-6343 Rotkreuz
Switzerland

Distribution

Roche Diagnostics GmbH
Sandhofer Straße 116
D-68305 Mannheim
Germany

Distribution in the US

Roche Diagnostics
9115 Hague Road
PO Box 50457
Indianapolis, IN 46250
USA

Trademarks

LIGHTCYCLER and LC are trademarks of Roche.

Other brands or product names are trademarks of their respective holders.

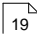
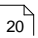
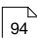
1.4 Glossary

EXOR	The database used by the LightCycler software. Multiple databases can be used on one instrument. This is to be avoided when using the LIMS software. see Other known problems ¹⁹
Experiment	A run on the instrument that results into various measurements (data points) being taken (Fluorescence and temperature data). These data are stored into the database and can be further analysed using various SW analyzing modules. The experiment data can be called via the LIMS interface.
IXO file	The format used by Exor to import/export objects like experiment results or other Database entities.
Macro	A set of parameters used to run an experiment. A macro is an Exor database object. The name of a valid macro must be passed to the LIMS interface to be able to start an experiment
MWP = Micro Well Plate	The container used to hold the samples to be analysed. The plates exist in two sizes, 96 (8x12) and 384 (16x24) wells. The samples are identified with coordinates : The 96-wells plate is numbered from A1 to H12; Rows A to H , columns 1 to 12. The 384-wells plate is numbered from A1 to P24; Rows A to P, columns 1 to 24.

2 LIMSCientLib Reference

This section of help provides reference information for the API elements provided by LIMSCientLib.

The LIMSCientLib reference information can be divided into following categories.

- Restrictions / Known problems  19
- Interfaces  20
- Enumerated Types  94

2.1 General information and debugging

2.1.1 General information

***** Conventions used in the manual *****

In some instances, the code on one line might be too long to fit properly on the page. To prevent the line from extending beyond pages margins, it has been split in the following way:

```
Public Function Login(PW As String, Host As String, ->
Port As String) As LIMSOperationResult
```

The -> sign indicates that the next line must be joined to recreate functional code, as in :

```
Public Function Login(PW As String, Host As String, Port As String) As LIMSOperationResult
```

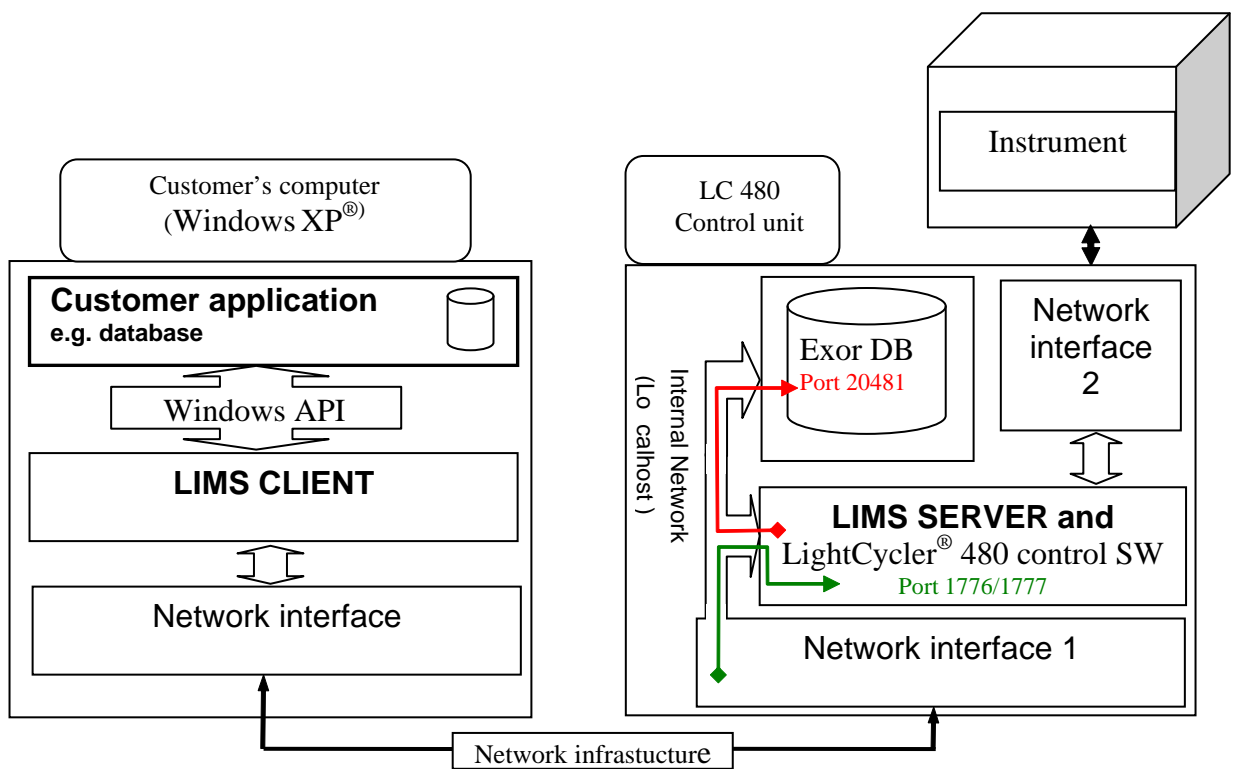
Important notes:

- Unless otherwise noted, all LIMS API calls are **synchronous**.
This means that the commands called will always wait until the command completes before control is returned to the calling program.
- The LIMS server expects to hear a **heartbeat** from a client every ten seconds.
For debugging purposes, you may want to set a break point and pause the client. To avoid getting logged off:
 1. Open the file C:\Program Files\Roche\LightCycler480\Bin\LIMS.stc
 2. Find and change the line:


```
<LIMS_SOCKET_SERVER_HEARTBEAT_TIMEOUT>10</LIMS_SOCKET_SERVER_HEARTBEAT_TIMEOUT>
```
 3. Replace ten with a larger value.
- The database being accessed by the LIMS server is the last one that was opened by the regular GUI software.
When using LIMS, it is a good practice to have only one operational database, otherwise there might be an uncertainty regarding the currently used database.
- If a firewall or other security device is present between the LIMS client and server, ports 1776 and 1777 must be opened to enable communication.

General architecture of the system:

Note that the Exor DB for the instrument might be installed at another location. e.g. Central server.



2. Connect from the application side.

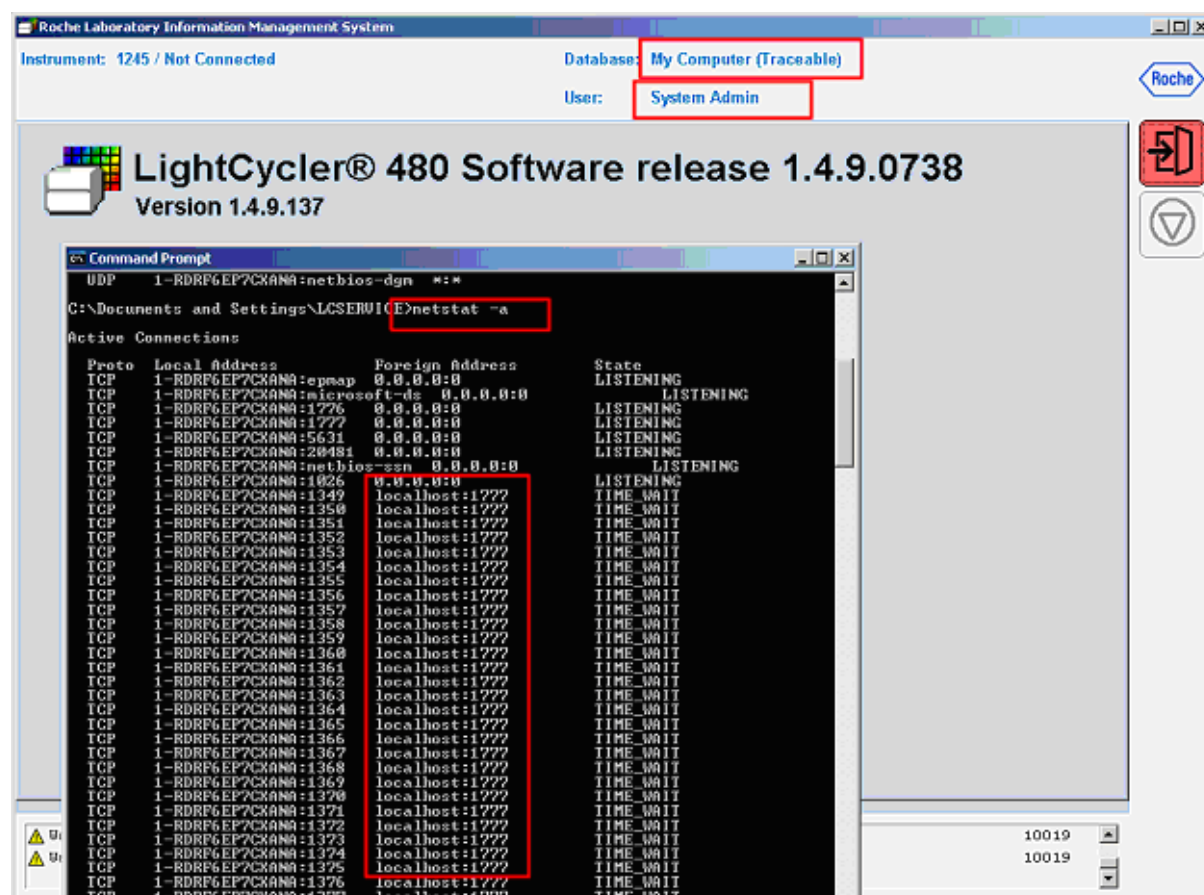
- The application should be on a separate computer running windows
- The LIMS client (Library) must be installed.

Once the application has successfully logged in and connected to the database, the following shows up on the LIMS server side:

- The database name. This is the last database that was successfully opened by the LightCycler application software.
- The user name that was used to connect.

If you retype the netstat -a command in the command window, you should probably see a lot of activity on the port 1777.

Additionally, if an active instrument has been successfully connected, it will show up on the LIMS Server screen with some status (here in the example, the default instrument was not able to connect).



2.1.2 Install LIMS server on another PC

To do read-only LIMS development (no instrument control) without impacting the instrument PC:

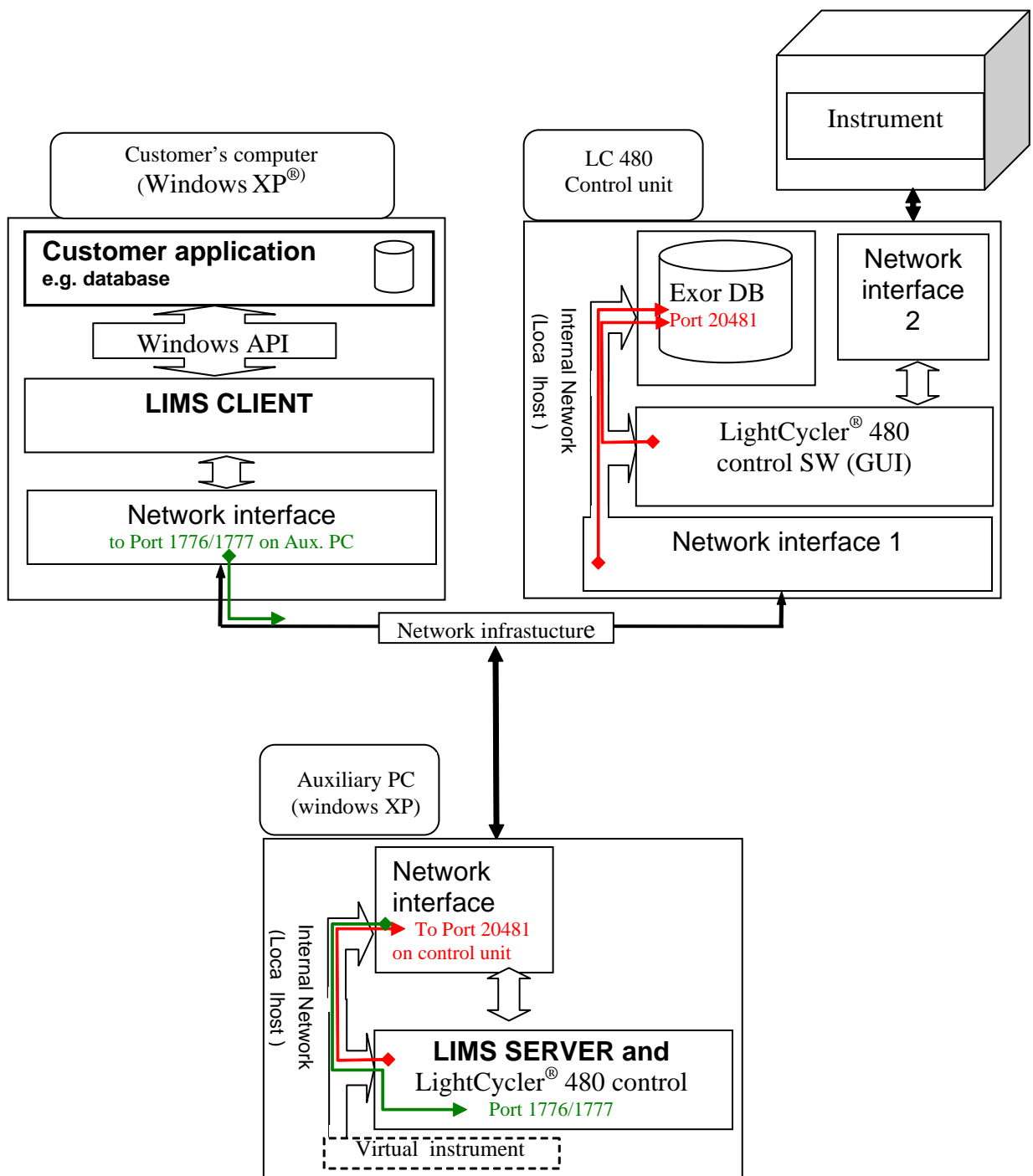
1. Install LC480 software and LIMS software on a second computer. (Named Auxiliary PC on the schematic)
2. On the Auxiliary PC start LSC480 software. On the login screen, hit the options button and configure a new database that points to Exor on the instrument PC (Control unit).
3. In LCS 480 software on the Auxiliary PC, set the default instrument to be a virtual instrument.
4. On the Auxiliary PC, log out of the LCS480 software and start the LIMS server
5. Connect a LIMS client to the LIMS server running on the Auxiliary PC

The LIMS client should be able to pull production data from the Exor on the instrument PC without affecting the operation of the instrument.

In other words, lab personnel should be able to run the instrument and informatics personnel should be able to retrieve data at the same time.

Of course, any instrument calls made by the LIMS client will fail, but queries for experiment data should work fine.

See schematic next page



2.1.3 Multiple LIMS users

Enabling multiple LIMS users on one server.

To enable multiple LIMS users on a single server:

1. Find the file

C:\Program Files\Roche\LightCycler480\Bin\LIMS.stc

2. In the file, find the line

```
<LIMS_SOCKET_SERVER_MAX_SESSIONS>1</LIMS_SOCKET_SERVER_MAX_SESSIONS>
```

3. Change the "1" to the number of user you want enable

```
<LIMS_SOCKET_SERVER_MAX_SESSIONS>4</LIMS_SOCKET_SERVER_MAX_SESSIONS>
```

Remember, only one user can access the instrument at a time.

2.2 Restrictions and known issues

2.2.1 Environment

The LIMS client must be installed on a computer with the following recommended characteristics:

Operating system:

Windows 2000 SP 4

Windows XP SP 2

2.2.2 Other known issues

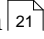
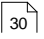
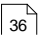
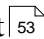
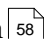
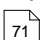

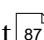

- The Database opened by the LIMS server by default is the last database opened by the regular GUI software.
Caution should be exercised if more than one database exist on the system.
For this reason it is **recommended** to install **only one database** when using the LIMS module,
- The Replicate Position property will be deprecated in the v 1.5. It is strongly recommended not to use it to ensure future compatibility.

Corrected issues in version 1.5

- The number of sockets is strictly enforced (see Multiple LIMS users¹⁸)
- Reservation of the instrument enforces exclusive access: If multiple sockets are enabled (see Multiple LIMS users¹⁸) the instrument can be reserved only by one socket at a time, even if the same user credentials are used.
When a connection has reserved the instrument, subsequent connections give only access to the database.
- Windows Language settings can be other than US. English.

2.3 Interfaces

This section contains reference information for the COM interfaces provided by LIMSCientLib. The following interfaces are used with LIMSCientLib.

- ILIMSConnection  21
- LIMSExperimentInfo  30
- LIMSIInstrument  36
- LIMSOperationResult  53
- LIMSSampleDefinition  58
- LIMSSampleInfo  71
- LIMSTQuery  80
- LIMSTQueryResult  87
- LIMSTQueryResultData  90

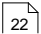
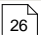
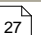
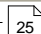
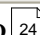
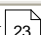
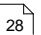
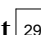
Requirements for all interfaces

Header: Declared in LIMSCientLib.h

Import Library: Use LIMSCientLib.lib

2.3.1 ILIMSConnection

Root interface for LIMS clients.

General methods and properties.		
property get	LoggedIn 	Checks to see if you are currently logged in to the LIMS server
property get / put	Host 	Gets the IP Address or machine name of the computer on which the LIMS server is running
property get / put	Port 	Gets the TCP port number on which the LIMS server is running (default = 1776)
property get	Instrument 	Returns an interface for controlling the instrument
property get	ExperimentInfo 	Returns an interface for accessing the Experiment API
property get	Query 	Returns an ILIMSQuery interface
method	Login 	Establishes a session with the LIMS server if not already logged in
method	Logout 	Terminates the session with the LIMS server

In order to use the LIMS interface, a **connection object** has to be created first. As an example, see the code fragments below:

[Visual Basic]

```
Public Sub CreateConnection()  
    If Not Assigned(gConnection) Then  
        Set gConnection = New LIMSClientLib.LIMSConnection  
        MsgBox ("Connection created")  
    Else  
        MsgBox ("Failure TO CREATE CONNECTION")  
    End If  
End Sub
```

[C#]

```
private bool CreateConnection()  
{  
    if (!Assigned("Object is already assigned"))  
    {  
        mConnection = new LIMSConnection();  
        VerifyMessage(login, "Connection created");  
        return true;  
    }  
    return false;  
}
```

2.3.1.1 LoggedIn Property

Checks to see if you are currently logged in to the LIMS server.

[C/C++]

```
HRESULT    get_LoggedIn( VARIANT_BOOL* Value );
```

[Visual Basic]

```
Public ReadOnly Property    LoggedIn As Boolean
```

Code Snippet

```
Public Function Logout() As LIMSOperationResult
On Error GoTo Error
    If gConnection.LoggedIn Then
        Set Logout = gConnection.Logout
        If Logout.Successful Then
            MsgBox ("Successfully Logged Out")
        Else
            MsgBox ("Failure to logout")
        End If
    End If
Error:
    Exit Function
End Function
```

[C#]

```
public ref bool    ILIMSConnection.LoggedIn { get; }
```

Parameters

Value

[out,retval]

2.3.1.2 Query Property

Returns an ILIMSQuery interface.

```
[C/C++]
HRESULT  get_Query( ILIMSQuery** Value );
```

```
[Visual Basic]
Public ReadOnly Property Query As Object
```

VB code snippet:

```
Public Function GetQueryResult()
    Dim QIndx, Indx As Integer
    Dim OpResult As LIMSOOperationResult
    Dim Qresult As LIMSQueryResult
    Dim Qdata As LIMSQueryResultData
    Dim LQuery As LIMSQuery

    If gConnection.LoggedIn = True Then
        gConnection.Query.ObjectType = "Experiment" 'set the filters
        SEE ALSO ObjectType
        gConnection.Query.Name = "*"
        gConnection.Query.Owner = ""
        Set OpResult = gConnection.Query.ExecuteQuery(Qresult)
        If OpResult.Successful = True Then
            QIndx = Qresult.Count 'Get the number of items returned
            If QIndx > 0 Then
                Main.ExpList.Clear ' Delete the list in main frame
                For Indx = 0 To QIndx - 1
                    'Read Individual items and Add them to the list
                    Set Qdata = Qresult.GetResultData(Indx)
                    Main.ExpList.AddItem (Qdata.Name), Indx
                Next Indx
            End If
        Else
            MsgBox "Exp Query unsuccessful .."
        End If
    End If
End Function
```

```
[C#]
public object  ILIMSConnection.Query { get; }
```

Parameters

Value

[out,retval]

2.3.1.3 ExperimentInfo Property

Returns an interface for accessing the Experiment API.

```
[C/C++]
HRESULT  get_ExperimentInfo( ILIMSExperimentInfo** Value );
```

```
[Visual Basic]
Public ReadOnly Property ExperimentInfo As Object
```

VB code snippet :

```
Public Function GetExperimentStatus(ExperimentName As String) ->
As LIMSEOperationResult
    Dim EStatus As String
    EStatus = ""
    If Assigned(gConnection) Then
Set GetExperimentStatus = gConnection.ExperimentInfo.GetStatus ->
(ExperimentName, EStatus)
        If GetExperimentStatus.Successful Then
            MsgBox ( "Status of experiment : " & EStatus)
        Else
            MsgBox ( "No Status found")
        End If
    Else
        MsgBox "Global LIMS Proxy object is not assigned"
    End If
End Function
```

```
[C#]
public object  ILIMSConnection.ExperimentInfo { get; }
```

Parameters

Value

[out,retval]

2.3.1.4 Instrument Property

Returns an interface for controlling the instrument.

```
[C/C++]  
HRESULT    get_Instrument( ILIMSIInstrument** Value );
```

```
[Visual Basic]  
Public ReadOnly Property Instrument As Object
```

VB code snippet :

```
Public Function ReserveInstrument() As LIMSOperationResult  
    If Assigned(gConnection) Then  
        Set ReserveInstrument = gConnection.Instrument.Reserve  
        If ReserveInstrument.Successful Then  
            MsgBox ("Instrument successfully reserved")  
        Else  
            MsgBox ("Failure")  
        End If  
    End If  
End Function
```

```
[C#]  
public object    ILIMSCConnection.Instrument { get; }
```

Parameters

Value

[out,retval]

2.3.1.5 Host Property

Gets the IP Address or machine name of the computer on which the LIMS server is running.

```
[C/C++]  
HRESULT  get_Host( BSTR* Value );  
HRESULT  put_Host( BSTR Value );
```

```
[Visual Basic]  
Public Overloads Property Host As String
```

VB code snippet :

```
Public Function Login(User As String, Password As String, Host As String, ->  
HostPort As String) As LIMSOperationResult  
    gConnection.Host = Host  
    gConnection.Port = HostPort  
    If Assigned(gConnection) Then  
        Set Login = gConnection.Login(User, Password)  
        If Login.Successful Then  
            MsgBox ("Logged in")  
        Else  
            MsgBox ("LoginFailure")  
        End If  
    Else  
        MsgBox "Global LIMS Proxy object is not assigned"  
    End If  
End Function
```

```
[C#]  
public ref string  ILIMSConnection.Host { get; set; }
```

Parameters

Value

[out,retval]

2.3.1.6 Port Property

Gets the TCP port number on which the LIMS server is listening (default = 1776).

[C/C++]

```
HRESULT get_Port( LONG* Value );  
HRESULT put_Port( LONG Value );
```

[Visual Basic]

```
Public Overloads Property Port As Long
```

VB code snippet: [see ILIMSConnection::Host](#)²⁶

[C#]

```
public ref int ILIMSConnection.Port { get; set; }
```

Parameters

Value

[out,retval]

2.3.1.7 Login Method

Establishes a login with the LIMS server if not already logged in.

```
[C/C++]
HRESULT Login(
    BSTR User,
    BSTR Password,
    ILIMSOperationResult** Value
);
```

```
[Visual Basic]
object.Login(
    ByVal User As String,
    ByVal Password As String,
) As Object
```

VB code snippet :

```
Public Function Login(User As String, Password As String, ->
Host As String, HostPort As String) As LIMSOOperationResult
    gConnection.Host = Host
    gConnection.Port = HostPort
    If Assigned(gConnection) Then
        Set Login = gConnection.Login(User, Password)
        If Login.Successful Then
            MsgBox ("Logged in")
        Else
            MsgBox ("LoginFailure")
        End If
    Else
        MsgBox "Global LIMS Proxy object is not assigned"
    End If
End Function
```

```
[C#]
void ILIMSConnection.Login(
    string User,
    string Password,
);
```

Parameters

User

[in]

Password

[in]

Value

[out,retval]

2.3.1.8 Logout Method

Terminates the login with the LIMS server.

```
[C/C++]
HRESULT Logout(
    LIMSOperationResult** Value
);
```

```
[Visual Basic]
object.Logout() As Object
```

VB code snippet :

```
Public Function Logout() As LIMSOperationResult
On Error GoTo Error
    If gConnection.LoggedIn Then
        Set Logout = gConnection.Logout
        If Logout.Successful Then
            MsgBox ("Successfully Logged Out")
        Else
            MsgBox ("Failure to logout")
        End If
    End If
Error:
    Exit Function
End Function
```

```
[C#]
object ILIMSConnection.Logout();
```

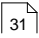
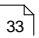

Parameters

Value

[out,retval]

2.3.2 ILIMSExperimentInfo

Interface providing information about the experiment.

General methods and properties.		
method	GetStatus 	Returns the current state of the specified experiment
method	GetCompletedExperimentSummary 	Generates and returns summary information for the specified experiment
method	ExportExperiment 	Exports an experiment to an IXO file

2.3.2.1 GetStatus Method

Returns the current state of the specified experiment.

```
[C/C++]
HRESULT GetStatus(
    BSTR ExperimentName,
    BSTR* Status,
    ILIMSOperationResult** Value
);
```

```
[Visual Basic]
object.GetStatus(
    ByVal ExperimentName as string,
    ByRef Status as string,
) as object
```

VB code snippet :

```
Public Function GetExperimentStatus(ExperimentName As String)->
As LIMSOOperationResult
    Dim EStatus As String
    EStatus = ""
    If Assigned(gConnection) Then
        Set GetExperimentStatus = gConnection.ExperimentInfo.->
GetStatus(ExperimentName, EStatus)
        If GetExperimentStatus.Successful Then
            MsgBox ( "Status of experiment : " & EStatus)
        Else
            MsgBox ( "No Status found")
        End If
    Else
        MsgBox "Global LIMS Proxy object is not assigned"
    End If
End Function
```

```
[C#]
void ILIMSExperimentInfo.GetStatus(
    string ExperimentName,
    ref string Status,
);
```

Parameters

ExperimentName
[in]

Status
[out]

Value
[out,retval]

Remarks

See table next page

The values for experiment status are as follows:

'No analyses':	experiment completed -- no analyses run
'Has analyses':	experiment completed -- at least one analysis run
'Not started':	experiment created, but not started in the instrument
'Running':	experiment running in the instrument
'Aborted':	experiment in the instrument was aborted
'Error':	experiment in the instrument had an error

2.3.2.2 GetCompletedExperimentSummary Method

Generates and returns summary information for the specified experiment.

[C/C++]

```
HRESULT GetCompletedExperimentSummary(  
    BSTR ExperimentName,  
    BSTR* Summary,  
    ILIMSOperationResult** Value  
);
```

[Visual Basic]

```
object.GetCompletedExperimentSummary(  
    ByVal ExperimentName As String,  
    ByRef Summary As String,  
) As Object
```

VB code snippet :

```
Public Function GetExperimentSummaryAsXML(ExpName)  
    Dim SummaryText As String  
    If Assigned(gConnection) Then  
        Set GetExperimentSummaryAsXML = gConnection.ExperimentInfo. ->  
GetCompletedExperimentSummary(ExpName, SummaryText)  
        If GetExperimentSummaryAsXML.Successful Then  
            Main.Text1.Text = SummaryText  
            ' To display: Use a large text container, the summary  
            ' can be > several tenth of MB  
            ' In VB a RichTextBox is suitable for such large objects.  
  
            MsgBox ("Size of summary =" & Str(Len(SummaryText)))  
            ' Length of the summary  
        Else  
            MsgBox ("Failure to acquire summary")  
        End If  
    Else  
        MsgBox " LIMS Proxy Not Assigned"  
    End If  
End Function
```

[C#]

```
void ILIMSExperimentInfo.GetCompletedExperimentSummary(  
    string ExperimentName,  
    ref string Summary,  
);
```

Parameters

ExperimentName

[in]

Summary

[out]

Value

[out,retval]See Notice next page

Notice

The experiment XML has several properties and four subsections, namely run, analysis, subsets, and samples:

```
<Experiment>
  <prop name="name">Standard-384-2</prop>
  more properties ...
  <run>
    run content ...
  </run>
  <analyses>
    analyses content ...
  </analyses>
  <Subsets>
    subset content ...
  <Subsets>
  <Samples>
    samples content ...
  </Samples>
</Experiment>
```

- Each section corresponds roughly to a section in the Experiment summary screen. The most commonly used section will be the analysis section, where analysis results can be obtained.
- For Absolute Quantification, Tm Calling, and Genotyping, each analysis section has a list of properties, followed by a list of samples. The analysis results for each sample are given in the list, and correspond to the columns in the result table of the analysis.
- For Relative Quantification, each analysis section has several properties and a list of pairings. For each pairing, the analysis results are listed as properties, and generally correspond to the columns on the result tab in the relative quantification analysis screen.

2.3.2.3 ExportExperiment Method

Exports an experiment to an IXO file.

[C/C++]

```
HRESULT    ExportExperiment(  
    BSTR          ExperimentName,  
    BSTR          Filename,  
    ILIMSOperationResult** Value  
);
```

[Visual Basic]

```
object.ExportExperiment(  
    ByVal ExperimentName As String,  
    ByVal Filename      As String,  
) As Object
```

[C#]

```
void    ILIMSExperimentInfo.ExportExperiment(  
    string ExperimentName,  
    string Filename,  
);
```

Parameters

ExperimentName

[in]

Filename

[in]

Value

[out,retval]

2.3.3 ILIMSIInstrument

Interface for controlling the instrument.

General methods and properties.		
method	Reserve ³⁷	Attempt to obtain exclusive control of the instrument
method	Unreserve ³⁸	Relinquish exclusive control of the instrument
method	Open ³⁹	Open the loading door on the instrument
method	Close ⁴¹	Close the loading door on the instrument
method	StartExperiment ⁴²	Run the specified experiment with the given parameters
method	GetStatus ⁴⁴	Returns the instrument status as a string
method	AbortExperiment ⁴⁶	Abort the currently running experiment
method	GetContainerBarcode ⁴⁷	Get the barcode of the tray currently loaded on the instrument
method	OpenAndWait ⁴⁸	open the loading door and wait for the status of the instrument
method	CloseAndWait ⁴⁹	close the loading door and wait for the status of the instrument

2.3.3.1 Reserve Method

Attempt to obtain exclusive control of the instrument

Before attempting to control an instrument, it should be reserved. Please note that exclusive control is only obtained if the relevant setup is correct.

See Multiple LIMS users¹⁸. Otherwise the control may be obtained by multiple clients.

[C/C++]

```
HRESULT Reserve(  
    ILIMSOperationResult** Value  
);
```

[Visual Basic]

```
object.Reserve() As Object
```

VB code snippet :

```
Public Function ReserveInstrument() As LIMSOperationResult  
    If Assigned(gConnection) Then  
        Set ReserveInstrument = gConnection.Instrument.Reserve  
        If ReserveInstrument.Successful Then  
            MsgBox "Instrument successfully reserved"  
        Else  
            ErNum = Str(ReserveInstrument.ErrorNumber)  
            ErMessage = ReserveInstrument.Message  
            MsgBox ( "Reserve failure info:" & ErNum & "/" & ErMessage)  
        End If  
    Else  
        MsgBox "Global LIMS Proxy object is not assigned"  
    End If  
End Function
```

[C#]

```
object ILIMSIInstrument.Reserve();
```

Parameters

Value

[out,retval]

2.3.3.2 Unreserve Method

Relinquish exclusive control of the instrument

[C/C++]

```
HRESULT Unreserve(  
    ILIMSOperationResult** Value  
);
```

[Visual Basic]

object.Unreserve() As Object

```
Public Function UnreserveInstrument() As LIMSOperationResult  
    If Assigned(gConnection) Then  
        Set UnreserveInstrument = gConnection.Instrument.Unreserve  
        If UnreserveInstrument.Successful Then  
            MsgBox ("Instrument successfully released")  
        Else  
            MsgBox ("Failure to unreserve")  
        End If  
    Else  
        MsgBox "Global LIMS Proxy object is not assigned"  
    End If  
End Function
```

[C#]

```
object ILIMSIInstrument.Unreserve();
```

Parameters

Value

[out,retval]

2.3.3.3 Open Method

Open the loading drawer on the instrument

Preliminary remarks.

As stated in the General information, the command is synchronous like most of the API methods.

When the command returns control a success condition will most probably occur, regardless of the mechanical movement being completed or not.

The LIMS API cannot assess the position of the plate drawer. The success condition refers only to the fact that the command was accepted and acknowledged by the system.

It is not possible to find directly from the success condition if the open movement was successfully completed or not, however, it is possible to find indirectly, by polling the instrument status after the open command completed "successfully".

If a mechanical failure occurred e.g. a mechanical obstacle prevented the movement of the drawer, the Instrument status will return an error condition

See GetStatus Method ⁴⁴

A successful Open should result in a '**Standby (no MWP)**' condition for the status testing.

Additionally, an attempt to open or close during the initialization phase of the analyzer will result into a timeout, but this will not be reported to the LIMS client, but rather as a "success".

For this reason, it is a good idea to test for the instrument standby condition before invoking the open function.

[C/C++]

```
HRESULT Open(  
    LIMSOperationResult** Value  
);
```

[Visual Basic]

```
object.Open() As Object
```

VB code snippet :

```
Public Function OpenInstrument() As LIMSOperationResult  
Dim TrayStatus As LIMSOperationResult  
If Assigned(gConnection) Then  
    Set OpenInstrument = gConnection.Instrument.Open  
    If OpenInstrument.Successful Then  
        MsgBox ( "Open command accepted"  
            ' Take some steps to check if the Instrument Status  
            ' is not 'Error'  
            ' See Function ILIMSIInstrument GetStatus44  
    Else  
        MsgBox ( "Open Command Failure"  
        ErNum = Str(OpenInstrument.ErrorNumber)  
        ErMessage = OpenInstrument.Message  
        MsgBox ^Box ( "Opening Error:" & ErNum & "/" & ErMessage)  
    End If  
Else  
    MsgBox "Global LIMS Proxy object is not assigned"  
End If  
End Function
```

[C#]

```
object ILIMSIInstrument.Open();
```

Parameters

Value

[out,retval]

2.3.3.4 Close Method

Close the loading drawer on the instrument

Preliminary remarks.

The same general remarks apply as for the Open function.

If a mechanical failure occurred e.g. a mechanical obstacle prevented the movement of the drawer, the Instrument status will return an error condition See GetStatus Method⁴⁴

A successful Close should result in a '**Standby (no MWP)**' or '**Standby (MWP loaded)**' condition for the status testing.

An attempt to open or close during the initialization phase of the analyzer, or if the drawer is already closed, will result into a timeout, but no error condition will be reported to the LIMS client in such cases.

[C/C++]

```
HRESULT Close(  
    ILIMSOperationResult** Value  
);
```

[Visual Basic]

object.Close() As Object

VB code snippet :

```
Public Function CloseInstrument() As LIMSOOperationResult  
    If Assigned(gConnection) Then  
        Set CloseInstrument = gConnection.Instrument.Close  
        If CloseInstrument.Successful Then  
            MsgBox ("Tray command successfully acknowledged")  
        Else  
            MsgBox ("Failure to close")  
        End If  
    Else  
        MsgBox "Global LIMS Proxy object is not assigned"  
    End If  
End Function
```

[C#]

```
object ILIMSIInstrument.Close();
```

Parameters

Value

[out,retval]

2.3.3.5 StartExperiment method

Run the specified experiment with the given parameters

See also the **SampleDefinition**⁵⁸ general Methods and properties.

- If sample definitions are to be used for the run, the collection of samples **must** be created (instantiated) **before** starting the experiment.
- A valid macro name must be passed to be able to start and experiment. The macro must exist and must belong to the logged in user, unless the user is administrator in the database.

[C/C++]

```
HRESULT StartExperiment(  
    BSTR ExperimentName,  
    BSTR ContainerBarCode,  
    BSTR MacroName,  
    ILIMSSampleDefinition* SampleDefinition,  
    ILIMSOperationResult** Value  
);
```

[Visual Basic]

```
object.StartExperiment(  
    ByVal ExperimentName As String,  
    ByVal ContainerBarCode As String,  
    ByVal MacroName As String,  
    ByRef SampleDefinition As ILIMSSampleDefinition,  
) As Object
```

VB code snippet :

```
Public Function StartExperiment(aExperimentName As String, aBarCodeContainer ->  
As String, aMacroName As String, aSampleDef As LIMSSampleDefinition) ->  
As LIMSOperationResult  
    Dim RetValue As Integer  
    If Assigned(gConnection) Then  
        Set StartExperiment = gConnection.Instrument.StartExperiment ->  
(aExperimentName, aBarCodeContainer, aMacroName, aSampleDef)  
        If StartExperiment.Successful Then  
            MsgBox ("Experiment started")  
        Else  
            ErNum = Str(StartExperiment.ErrorNumber)  
            ErMessage = StartExperiment.Message  
            RetValue = MsgBox ("Error:" & ErNum & "/" & ErMessage,vbOKOnly &  
vbInformation,"Experiment Start Failed")  
        End If  
    Else  
        MsgBox "Global LIMS Proxy object is not assigned"  
    End If  
End Function
```

[C#]

```
void ILIMSIInstrument.StartExperiment(  
    string ExperimentName,  
    string ContainerBarCode,  
    string MacroName,  
    ref ILIMSSampleDefinition SampleDefinition,  
);
```

Parameters

ExperimentName

[in]

ContainerBarCode

[in]

MacroName

[in]

SampleDefinition

[in]

Value

[out,retval]

2.3.3.6 GetStatus Method

Returns the instrument status as a string

Note : This command is also useful to find out if a **close** or **open** command was successfully completed. See Open Method^[39] and Close Method^[41]

[C/C++]

```
HRESULT GetStatus(  
    BSTR* Status,  
    ILIMSOperationResult** Value  
);
```

[Visual Basic]

```
object.GetStatus(  
    ByRef Status As String,  
) As Object
```

VB code snippet :

```
Public Function GetInstrumentStatus() As LIMSOperationResult  
    Dim IStatus As String  
    IStatus = ""  
    If Assigned(gConnection) Then  
        Set GetInstrumentStatus = gConnection.Instrument.GetStatus(IStatus)  
        If GetInstrumentStatus.Successful Then  
            Select Case IStatus  
                Case "Error"  
                    ' Action for the ERROR status e.g. Set a Global Flag....  
                Case "Initializing"  
                    ' Action for the INITIALIZING status  
                Case Else  
                    ' Action for the OTHER status  
            End Select  
            MsgBox ("Instrument Status: " & IStatus)  
        Else  
            MsgBox ("Failure to get instrument STATUS")  
        End If  
    Else  
        MsgBox "Global LIMS Proxy object is not assigned"  
    End If  
End Function
```

[C#]

```
void ILIMSIInstrument.GetStatus(  
    ref string Status,  
);
```

Parameters

Status

[out]

Value

[out,retval]

Remarks See table next page

The values for instrument status are as follows:

'Running':	The instrument is running an experiment
'Service':	The instrument is being serviced
'Error':	The instrument has an error
'Initializing':	The instrument is initializing
'Standby (MWP loaded)':	The instrument is ready with a plate loaded
'Standby (no MWP)':	The instrument is on standby with no plate loaded

2.3.3.7 AbortExperiment Method

Abort the currently running experiment
[C/C++]

```
HRESULT AbortExperiment(  
    ILIMSOperationResult** Value  
);
```

[Visual Basic]

object.AbortExperiment() As Object

VB code snippet :

```
Public Function AbortExperiment() As LIMSOperationResult  
Dim RetValue As Integer  
If Assigned(gConnection) Then  
    Set AbortExperiment = gConnection.Instrument.AbortExperiment  
    If AbortExperiment.Successful Then  
        MsgBox ("Abort experiment sucessful")  
        Main.CdeAbort.BackColor = &H8000000F  
    Else  
        ErNum = Str(AbortExperiment.ErrorNumber)  
        ErMessage = AbortExperiment.Message  
        RetValue = MsgBox (" Error:" & ErNum & "/" & ErMessage, vbOKOnly &  
            vbInformation , "Abort experiment Failed")  
    End If  
End If
```

[C#]

```
object ILIMSIInstrument.AbortExperiment();
```

Parameters

Value

[out,retval]

2.3.3.8 GetContainerBarcode Method

Get the barcode of the tray currently loaded on the instrument

[C/C++]

```
HRESULT GetContainerBarcode(  
    BSTR* Barcode,  
    ILIMSOperationResult** Value  
);
```

[Visual Basic]

```
object.GetContainerBarcode(  
    ByRef Barcode As String,  
) As Object
```

VB code snippet :

```
Public Function GetContainerBarcode() As ILIMSOperationResult  
    Dim CBarcode As String  
    If Assigned(gConnection) Then  
        Set GetContainerBarcode = gConnection.Instrument.->  
GetContainerBarcode(CBarcode)  
        If GetContainerBarcode.Successful Then  
            MsgBox ("Container Barcode: " & CBarcode)  
  
        Else  
            MsgBox ("Failure no BC found")  
        End If  
    Else  
        MsgBox "Global LIMS Proxy object is not assigned"  
    End If  
End Function
```

[C#]

```
void ILIMSIInstrument.GetContainerBarcode(  
    ref string Barcode,  
);
```

Parameters

Barcode

[out]

Value

[out,retval]

2.3.3.9 OpenAndWait Method

open the loading door and wait for the status of the instrument

Preliminary remark.

When the normal `open`³⁹ command returns control to the API a success condition will most probably occur regardless of the mechanical movement being completed or not. The LIMS API cannot assess the position of the plate drawer.

This method provides an alternative to `open + get status`.

As a successful Open should result in a '**Standby (no MWP)**' condition, this is normally the parameter that should be passed as the **StatusToWaitFor**.

[C/C++]

```
HRESULT OpenAndWait(  
    LONG TimeoutSeconds,  
    BSTR StatusToWaitFor,  
    ILIMSOperationResult** Value  
);
```

[Visual Basic]

```
object.OpenAndWait(  
    ByVal TimeoutSeconds as long,  
    ByVal StatusToWaitFor as string,  
) as object
```

[C#]

```
void ILIMSIInstrument.OpenAndWait(  
    int TimeoutSeconds,  
    string StatusToWaitFor,  
);
```

Parameters

TimeoutSeconds

[in]

StatusToWaitFor , normally '**Standby (no MWP)**'

[in]

Value

[out,retval]

2.3.3.10 CloseAndWait Method

close the loading door and wait for the status of the instrument

Preliminary remark.

When the normal `close`⁴¹ command returns control to the API a success condition will most probably occur regardless of the mechanical movement being completed or not. The LIMS API cannot assess the position of the plate drawer.

This method provides an alternative to `close + get status`.

As a successful Close should result in a '**Standby (no MWP)**' or '**Standby (MWP loaded)**' condition.

Unless the program has also full control on the loading of a plate, e.g. robotic arm, it is difficult to use this function as the correct expected condition is not completely predictable.

[C/C++]

```
HRESULT CloseAndWait(  
    LONG TimeoutSeconds,  
    BSTR StatusToWaitFor,  
    ILIMSOperationResult** Value  
);
```

[Visual Basic]

```
object.CloseAndWait(  
    ByVal TimeoutSeconds As long,  
    ByVal StatusToWaitFor As String,  
) As Object
```

[C#]

```
void ILIMSIInstrument.CloseAndWait(  
    int TimeoutSeconds,  
    string StatusToWaitFor,  
);
```

Parameters

TimeoutSeconds

[in]

StatusToWaitFor

[in]

Value

[out,retval]

2.3.3.11 GetContainerSensor

Get the current setting of the container (plate) sensor

Preliminary remarks.

The plate sensor detects the difference between MWP with either 96 or 384 wells. The sensor can function only with white plates.

When clear plates (transparent plastic) are used, the sensor must be turned OFF.

For the manual use of the instruments, the Administrator can setup general parameters in one of three way:

1. Always ON (White plates only)
2. Always OFF (Clear plates)
3. Mixed (The user decides at the start of each run)

Independent of those settings, the LIMS interface can turn the sensor ON or OFF at any time. (See also SetContainerSensor⁵²)

Values returned by the method:
can be "ON"; "OFF"; "failed"

[C/C++]

```
HRESULT GetContainerSensor(  
    BSTR* Sensor,  
    ILIMSOperationResult** Value  
);
```

[Visual Basic]

```
object.GetContainerSensor(  
    ByRef Sensor as string,  
) as object
```

VB code snippet :

```
Public Function GetContSens_value() As ILIMSOperationResult  
    Dim ContSensVal As String  
    If Assigned(gConnection) Then  
        Set GetContSens_value = ->  
            gConnection.Instrument.GetContainerSensor(ContSensVal)  
        Main.PlateType.Text = ContSensVal ' display result in a text box  
    Else  
        MsgBox "Global LIMS Proxy object is not assigned"  
    End If  
End Function
```

[C#]

```
void ILIMSIInstrument.GetContainerSensor(  
    ref string Sensor,  
);
```

Parameters

Sensor
[out]

[out,retval]

2.3.3.12 SetContainerSensor

Set the Container (plate) Sensor to the specified value

Value can be "ON"; "OFF"

```
[C/C++]
HRESULT SetContainerSensor(
    BSTR Sensor,
    ILIMSOperationResult** Value
);
```

[Visual Basic]

```
object.SetContainerSensor(
    ByVal Sensor As String,
) As Object
```

VB code snippet :

```
Public Function SetContSens_value(SensVal As String)
    ' SensVal can be 'ON' or 'OFF'
    If Assigned(gConnection) Then
        gConnection.Instrument.SetContainerSensor (SensVal)
    Else
        MsgBox "Global LIMS Proxy object is not assigned"
    End If
End Function
```

[C#]

```
void ILIMSIInstrument.SetContainerSensor(
    string Sensor,
);
```

Parameters

Sensor

[in]

Value

[out,retval]

2.3.4 ILIMSOOperationResult

Dispatch interface for LIMSOOperationResult Object

General methods and properties.		
property get	Successful ⁵⁴	True if the operation was successful, false otherwise
property get	ServerError ⁵⁵	True if error occurred on LIMS server rather than client
property get	Message ⁵⁵	Technical error message
property get	ErrorNumber ⁵⁶	Error number
property get	DateTime ⁵⁶	Date and time that the error occurred
property get	UserMessage ⁵⁷	End user error message

2.3.4.1 Successful Property

True if the operation was successful, false otherwise

```
[C/C++]
HRESULT get_Successful( VARIANT_BOOL* Value );
```

```
[Visual Basic]
Public ReadOnly Property Successful As Boolean
```

VB code snippet :

```
.....
Set ReserveInstrument = gConnection.Instrument.Reserve
If ReserveInstrument.Successful Then
    MsgBox "Instrument successfully reserved"
Else
    If ReserveInstrument.ServerError then          'Detect SERVER errors
        Stamp = ReserveInstrument.DateTime         'Server error Time Stamp
        ErNum = Str(ReserveInstrument.ErrorNumber)'Fetch Error number
        ErMessage = ReserveInstrument.Message      'Fetch Message
        RetValue = MsgBox ( "Err Number : " & ErNum & "/" & _
            ErMessage,vbOKOnly & vbInformation, "Reserve failure info at " &Tstamp )
    End If
End If
.....
```

```
[C#]
public ref bool ILIMSOperationResult.Successful { get; }
```

Parameters

Value

[out,retval]

2.3.4.2 ServerError Property

True if error occurred on LIMS server rather than client

```
[C/C++]  
HRESULT    get_ServerError( VARIANT_BOOL* Value );
```

```
[Visual Basic]  
Public ReadOnly Property ServerError As Boolean
```

VB code snippet : see Successful Property⁵⁴

```
[C#]  
public ref bool    ILIMSOperationResult.ServerError { get; }
```

Parameters

Value
[out,retval]

2.3.4.3 Message Property

Technical error message

```
[C/C++]  
HRESULT    get_Message( BSTR* Value );
```

```
[Visual Basic]  
Public ReadOnly Property Message As String
```

VB code snippet : see Successful Property⁵⁴

```
[C#]  
public ref string    ILIMSOperationResult.Message { get; }
```

Parameters

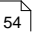
Value
[out,retval]

2.3.4.4 ErrorNumber Property

Error number

```
[C/C++]
HRESULT  get_ErrorNumber( LONG* Value );
```

```
[Visual Basic]
Public ReadOnly Property  ErrorNumber As Long
```

VB code snippet : see Successful Property 

```
[C#]
public  ref int  ILIMSOperationResult.ErrorNumber { get; }
```

Parameters

Value

[out,retval]

2.3.4.5 DateTime Property

Date and time when the error occurred (according to the LIMS CLIENT clock)

```
[C/C++]
HRESULT  get_DateTime( DATE* Value );
```

```
[Visual Basic]
Public ReadOnly Property  DateTime As Date
```

VB code snippet : see Successful Property 

```
[C#]
public  ref System.DateTime  ILIMSOperationResult.DateTime { get; }
```

Parameters

Value

[out,retval]

2.3.4.6 UserMessage Property

End user error message

Note : This property does not yield any content in current version.

```
[C/C++]  
HRESULT    get_UserMessage( BSTR* Value );
```

```
[Visual Basic]  
Public ReadOnly Property UserMessage As String
```

```
[C#]  
public ref string    ILIMSOperationResult.UserMessage { get; }
```

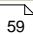
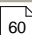
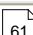
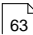
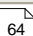
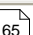
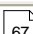
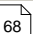
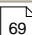
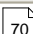
Parameters

Value

[out,retval]

2.3.5 ILIMSSampleDefinition

Represents a collection of samples

General methods and properties.		
property get	SampleCount 	returns the number of samples in the collection
property get	ExpDataCount 	Returns the number of experiment data strings
method	AddSample 	Adds a sample to the collection
method	Clear 	Removes all samples from the collection
method	DeleteSample 	Delete the sample at the given index position
method	GetSample 	Returns the sample at the given index position
method	AddExpData 	Add an experiment data string
method	GetExpData 	Return the experiment data string at the given index
method	DeleteExpData 	Delete the experiment data string at the given index
method	ClearExpData 	Delete all experiment data strings

A Sample collection (ILIMSSampledefinition) object must be created before the collection of samples can be manipulated by the program:

[Visual Basic]

VB code snippet :

```
set gSampleDef = New LIMSSampleDefinition
```

2.3.5.1 SampleCount Property

returns the number of samples in the collection

[C/C++]

```
HRESULT get_SampleCount( LONG* Value );
```

[Visual Basic]

```
Public ReadOnly Property SampleCount As Long
```

```
' The sample collection gSampleDef must exist prior to calling this routine.  
' e.g. as a Global variable. see ILIMSSampleDefinition58
```

```
Public Sub GetSampleCount()  
    Dim SampCount As Long  
    On Error GoTo Error  
    If (Assigned(gConnection) And gConnection.LoggedIn = True) Then  
        SampCount = gSampleDef.SampleCount  
        RetValue = MsgBox(Str(SampCount) &_  
            " Samples found", vbOKOnly, " Samples in collection")  
    End If  
Exit Sub
```

[C#]

```
public ref int ILIMSSampleDefinition.SampleCount { get; }
```

Parameters

Value

[out,retval]

2.3.5.2 ExpDataCount Property

Returns the number of experiment data strings

```
[C/C++]  
HRESULT  get_ExpDataCount( LONG* Value );
```

```
[Visual Basic]  
Public ReadOnly Property  ExpDataCount As Long
```

```
[C#]  
public  ref int  ILIMSSampleDefinition.ExpDataCount { get; }
```

Parameters

Value

[out,retval]

2.3.5.3 AddSample Method

Adds a sample to the collection

Notice: The Replicate Position property will be deprecated in the v 1.5.
It is strongly recommended not to use it to ensure future compatibility

```
[C/C++]
HRESULT AddSample(
    ILIMSSampleInfo* Sample
);
```

```
[Visual Basic]
object.AddSample(
    ByRef Sample As ILIMSSampleInfo
```

VB code snippet :

```
Public Sub AddSample()
' Creates the sample infos and populates them with the sample data.
' Then adds each of them to the sample definition (this object must have been
' instantiated before See ILIMSSampleDefinition58)

    Dim gSampleInfo As ILIMSSampleInfo
    Dim SampCount As Long
    Dim SamplePos As String
    gSampleDef.Clear ' First clear the whole collection

    'The sample set must be first instantiated
    Set gSampleInfo = New LIMSSampleInfo
    gSampleInfo.Position = "A1"
    gSampleInfo.ID = "FirstID"
    gSampleInfo.Name = "Sample Name XYZ"
    gSampleInfo.Comment = "Dummy Sample"

    gSampleDef.AddSample gSampleInfo ' Add the set to the collection.

    ' Set a sample info for a replicate
    Set gSampleInfo = New LIMSSampleInfo

    gSampleInfo.Position = "A3"
    gSampleInfo.ID = "" ' Name, ID, Comment MUST be empty for repl. pos.
    gSampleInfo.Name = ""
    gSampleInfo.Comment = ""
    gSampleInfo.ReplicatePosition = "A1"

    gSampleDef.AddSample gSampleInfo

End Sub
)
```

[C#]

```
void    ILIMSSampleDefinition.AddSample(  
    ref ILIMSSampleInfo Sample  
);
```

Parameters

Sample

[in]

2.3.5.4 Clear Method

Removes all samples from the collection.

[C/C++]

```
HRESULT Clear();
```

[Visual Basic]

```
object.Clear()
```

VB code snippet : see AddSample Method 

[C#]

```
void ILIMSSampleDefinition.Clear();
```

2.3.5.5 DeleteSample Method

Delete the sample at the given index position

[C/C++]

```
HRESULT DeleteSample(  
    LONG Index  
);
```

[Visual Basic]

```
object.DeleteSample(  
    ByVal Index As Long  
)
```

[C#]

```
void ILIMSSampleDefinition.DeleteSample(  
    int Index  
);
```

Parameters

Index

[in]

2.3.5.6 GetSample Method

Returns the sample at the given index position

[C/C++]

```
HRESULT GetSample(  
    LONG Index,  
    ILIMSSampleInfo** Value  
);
```

[Visual Basic]

```
object.GetSample(  
    ByVal Index As Long,  
) As Object
```

```
Public Sub SampFback(Scount As Integer) ' Scount = Number of samples found  
                                         ' in the collection (See SampleCount Property) 59  
    Dim Disp As String  
    Dim i, j, X As Integer  
    ' SampFB is a RichTextBox control to display the list with tabs  
    ' Samp_Feedback is the v.b. form containing the above RichTextBox  
  
    ScaleMode = vbPixels  
    Disp = ""  
    Samp_Feedback.Show  
    SampFB.Text = Disp  
  
    SampFB.SelStart = 0  
    SampFB.SelLength = Len(Samp_Feedback.SampFB.Text)  
    SampFB.SelTabCount = 4  
    SampFB.SelTabs(0) = 50  
    SampFB.SelTabs(1) = 200  
    SampFB.SelTabs(2) = 300  
    SampFB.SelTabs(3) = 600  
  
    Disp = "Position" & vbTab & "ID" & vbTab & _  
           "Name" & vbTab & "Comment" & vbCrLf & vbCrLf  
    SampFB.Text = Disp  
  
    For i = 0 To Scount - 1  
        Set gSampleInfo = gSampleDef.GetSample(i)  
        Disp = Disp & gSampleInfo.Position & vbTab & gSampleInfo.ID & _  
               vbTab & gSampleInfo.Name & vbTab & gSampleInfo.Comment & vbCrLf  
        SampFB.Text = Disp  
    Next i  
Exit Sub  
Error:  
RetVal = MsgBox("Error#" & Str(Err.Number) & _  
vbCrLf & Err.Description, vbOKOnly, "Error sample data read")  
End Sub
```

```
[C#]  
void    ILIMSSampleDefinition.GetSample(  
    int    Index,  
);
```

Parameters

Index

[in]

Value

[out,retval]

2.3.5.7 AddExpData Method

Add an experiment data string

[C/C++]

```
HRESULT AddExpData(  
    BSTR ExpData  
);
```

[Visual Basic]

```
object.AddExpData(  
    ByVal ExpData As String  
)
```

[C#]

```
void ILIMSSampleDefinition.AddExpData(  
    string ExpData  
);
```

Parameters

ExpData

Remarks

In LightCycler 480 version 1.5, the experiment data object is used only to set the external experiment name for a relative quantification or PCR Endpoint Genotyping experiments.

Relative Quantification uses the external experiment name as well as the subset and program. The subset name should be the same name that is displayed in subset editor.

The program number is a zero-based integer referring to programs as listed in the experiment. The XML string that should be passed to AddExpData is

```
<RelQuant>  
  <ReferenceExperimentName>experiment name goes here</ReferenceExperimentName>  
  <ReferenceExperimentSubsetName>subset name goes  
here</ReferenceExperimentSubsetName>  
  <ReferenceExperimentProgramNumber>nprogram number goes here... first program is  
0</ReferenceExperimentProgramNumber>  
</RelQuant>
```

Endpoint Genotyping requires only the experiment name. The XML string that should be passed to AddExpData is

```
<EndptGeno>  
  <ReferenceExperimentName>name goes here</ReferenceExperimentName>  
</EndptGeno>
```

Note that the XML values (i.e. "name goes here") must be replaced by actual values..

2.3.5.8 GetExpData Method

Return the experiment data string at the given index

[C/C++]

```
HRESULT GetExpData(  
    LONG Index,  
    BSTR* Value  
);
```

[Visual Basic]

```
object.GetExpData(  
    ByVal Index As Long,  
) As String
```

[C#]

```
void ILIMSSampleDefinition.GetExpData(  
    int Index,  
);
```

Parameters

Index

[in]

Value

[out,retval]

2.3.5.9 DeleteExpData Method

Delete the experiment data string at the given index

[C/C++]

```
HRESULT DeleteExpData(  
    LONG Index  
);
```

[Visual Basic]

```
object.DeleteExpData(  
    ByVal Index As Long  
)
```

[C#]

```
void ILIMSSampleDefinition.DeleteExpData(  
    int Index  
);
```

Parameters

Index

[in]

2.3.5.10 ClearExpData Method

Delete all experiment data strings

[C/C++]

```
HRESULT ClearExpData();
```

[Visual Basic]

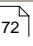
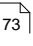
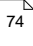
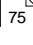
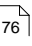
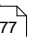
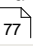
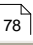
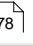
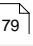
```
object.ClearExpData()
```

[C#]

```
void ILIMSSampleDefinition.ClearExpData();
```

2.3.6 ILIMSSampleInfo

Represents a single sample in a SampleDefinition collection

General methods and properties.		
property get / put	Position 	Get the sample position
property get / put	Name 	Get the sample name
property get / put	ID 	Get the sample ID
property get / put	Comment 	Get the sample comment
property get / put	ReplicatePosition 	Get the address of the master sample (i.e. A1 or C3) for replicates
property get	ExpDataCount 	(for future use) Returns the number of experiment data strings for this sample
method	AddExpData 	(for future use) Add an experiment data string for this sample
method	GetExpData 	(for future use) Return the experiment data string at the given index
method	DeleteExpData 	(for future use) Delete the experiment data string at the given index
method	ClearExpData 	(for future use) Delete all experiment data strings

2.3.6.1 Position Property

Get/Put the sample position.

```
[C/C++]  
HRESULT  get_Position( BSTR* Value );  
HRESULT  put_Position( BSTR Value );
```

```
[Visual Basic]  
Public Overloads Property Position As String
```

VB code snippet : see AddSample Method 

```
[C#]  
public ref string  ILIMSSampleInfo.Position { get; set; }
```

Parameters

Value

[out,retval]

2.3.6.2 Name Property

Get the sample name or assign it to the Name property.

Note : when creating a sample position, the sample name **cannot** be blank ("", or empty string).

```
[C/C++]
HRESULT  get_Name( BSTR* Value );
HRESULT  put_Name( BSTR Value );
```

```
[Visual Basic]
Public Overloads Property Name As String
```

VB code snippet : see AddSample Method 

```
[C#]
public ref string  ILIMSSampleInfo.Name { get; set; }
```

Parameters

Value

[out,retval]

2.3.6.3 ID Property

Get the sample ID

```
[C/C++]  
HRESULT  get_ID( BSTR* Value );  
HRESULT  put_ID( BSTR Value );
```

```
[Visual Basic]  
Public Overloads Property ID As String
```

VB code snippet : see AddSample Method 

```
[C#]  
public ref string  ILIMSSampleInfo.ID { get; set; }
```

Parameters

Value

[out,retval]

2.3.6.4 Comment Property

Get the sample comment

```
[C/C++]  
HRESULT get_Comment( BSTR* Value );  
HRESULT put_Comment( BSTR Value );
```

```
[Visual Basic]  
Public Overloads Property Comment As String
```

VB code snippet : see AddSample Method⁶¹

```
[C#]  
public ref string ILIMSSampleInfo.Comment { get; set; }
```

Parameters

Value

[out,retval]

2.3.6.5 ReplicatePosition Property *****

Get the address of the master sample (i.e. A1 or C3) for replicates

Notice: The Replicate Position property will be deprecated in the v 1.5.
It is strongly recommended not to use it to ensure future compatibility

```
[C/C++]  
HRESULT    get_ReplicatePosition( BSTR* Value );  
HRESULT    put_ReplicatePosition( BSTR Value );
```

```
[Visual Basic]  
Public Overloads Property ReplicatePosition As String
```

VB code snippet : see AddSample Method 

```
[C#]  
public ref string    ILIMSSampleInfo.ReplicatePosition { get; set; }
```

Parameters

Value

[out,retval]

2.3.6.6 ExpDataCount Property

Returns the number of experiment data strings for this sample

```
[C/C++]
HRESULT    get_ExpDataCount( LONG* Value );
```

```
[Visual Basic]
Public ReadOnly Property ExpDataCount As Long
```

```
[C#]
public ref int    ILIMSSampleInfo.ExpDataCount { get; }
```

Parameters

Value
[out,retval]

2.3.6.7 AddExpData Method

Add an experiment data string for this sample

```
[C/C++]
HRESULT    AddExpData(
    BSTR ExpData
);
```

```
[Visual Basic]
object.AddExpData(
    ByVal ExpData As String
)
```

```
[C#]
void    ILIMSSampleInfo.AddExpData(
    string ExpData
);
```

Parameters

ExpData
[in]

2.3.6.8 GetExpData Method

Return the experiment data string at the given index.

```
[C/C++]
HRESULT GetExpData(
    LONG Index,
    BSTR* Value
);
```

```
[Visual Basic]
object.GetExpData(
    ByVal Index As Long,
) As String
```

```
[C#]
void ILIMSSampleInfo.GetExpData(
    int Index,
);
```

Parameters

Index

[in]

Value

[out,retval]

2.3.6.9 DeleteExpData Method

Delete the experiment data string at the given index.

```
[C/C++]
HRESULT DeleteExpData(
    LONG Index
);
```

```
[Visual Basic]
object.DeleteExpData(
    ByVal Index As Long
)
```

```
[C#]
void ILIMSSampleInfo.DeleteExpData(
    int Index
);
```

Parameters

Index

[in]

2.3.6.10 ClearExpData Method

Delete all experiment data strings.

[C/C++]

```
HRESULT ClearExpData();
```

[Visual Basic]

```
object.ClearExpData()
```

[C#]

```
void ILIMSSampleInfo.ClearExpData();
```

2.3.7 ILIMSQuery

General methods and properties.		
property get / put	ObjectType ⁸¹	The type of object to return, or empty string for all object types.
property get / put	Name ⁸²	Filter for object name
property get / put	Owner ⁸²	Filter for owner of object
property get / put	FromDate ⁸³	Begin date search range for object creation or modification date
property get / put	ToDate ⁸⁴	End date search range for object creation or modification date
property get / put	QueryDateType ⁸⁵	Type of date search, either a creation date, modification date, or both. See the enumerated type for LIMSQueryDateType.
method	ExecuteQuery ⁸⁶	Execute the query

See Also [ILIMSQueryResultData::Path⁹³](#)

2.3.7.1 ObjectType Property

The type of object to return, or empty string for all object types.

[C/C++]

```
HRESULT get_ObjectType( BSTR* Value );
HRESULT put_ObjectType( BSTR Value );
```

[Visual Basic]

```
Public Overloads Property ObjectType As String
```

VB code snippet :

```
Public Function GetQueryResult()
    Dim QIndx, Indx As Integer
    Dim OpResult As LIMSOperationResult
    Dim Qresult As LIMSQueryResult
    Dim Qdata As LIMSQueryResultData
    Dim LQuery As LIMSQuery

    If gConnection.LoggedIn = True Then
        gConnection.Query.ObjectType = "Experiment" 'set the filters
                                                    SEE Remarks below

        gConnection.Query.Name = ""
        gConnection.Query.Owner = ""
        Set OpResult = gConnection.Query.ExecuteQuery(Qresult)
        If OpResult.Successful = True Then
            QIndx = Qresult.Count 'Get the number of items returned
            If QIndx > 0 Then
                Main.ExpList.Clear ' Delete the list in main frame
                For Indx = 0 To QIndx - 1
                    'Read Individual items and Add them to the list
                    Set Qdata = Qresult.GetResultData(Indx)
                    Main.ExpList.AddItem (Qdata.Name), Indx
                Next Indx
            End If
        Else
            MsgBox "Exp Query unsuccessful .."
        End If
    End If
End Function
```

[C#]

```
public ref string ILIMSQuery.ObjectType { get; set; }
```

Parameters

Value

[out,retval]

Remarks

Valid object types are as follows: Macro, Experiment, ColorComp, StdCurve, Template

See Also [ILIMSQueryResultData::Path⁹³](#)

2.3.7.2 Name Property

Filter for object name

```
[C/C++]
HRESULT  get_Name( BSTR* Value );
HRESULT  put_Name( BSTR Value );
```

```
[Visual Basic]
Public Overloads Property Name As String
```

VB code snippet	See	ILIMSTQuery::ObjectType ⁸¹
	See Also	ILIMSTQueryResultData::Path ⁹³

```
[C#]
public ref string  ILIMSTQuery.Name { get; set; }
```

Parameters

Value
[out,retval]

2.3.7.3 Owner Property

Filter for owner of object

```
[C/C++]
HRESULT  get_Owner( BSTR* Value );
HRESULT  put_Owner( BSTR Value );
```

```
[Visual Basic]
Public Overloads Property Owner As String
```

VB code snippet	See	ILIMSTQuery::ObjectType ⁸¹
-----------------	-----	---

```
[C#]
public ref string  ILIMSTQuery.Owner { get; set; }
```

Parameters

Value
[out,retval]

2.3.7.4 FromDate property

Begin date search range for object creation or modification date

```
[C/C++]
HRESULT get_FromDate( DATE* Value );
HRESULT put_FromDate( DATE Value );
```

```
[Visual Basic]
Public Overloads Property FromDate As Date
```

VB code snippet

```
...
Dim StartDate, EndDate As Date
StartDate = "6/9/2005"
EndDate = "6/10/2005"
...
gConnection.Query.Name = "*"
gConnection.Query.Owner = ""
gConnection.Query.QueryDateType = qdtCreationDateQuery ' Created between
gConnection.Query.FromDate = StartDate                  ' First modified
gConnection.Query.ToDate = EndDate                      ' Last modified
    Set OpResult = gConnection.Query.ExecuteQuery(Qresult)
...
gConnection.Query.QueryDateType = qdtAllDateQuery ' Reset the date if you
                                                    need to run further
                                                    queries w/o dates
```

See also [ILIMSQuery::ObjectType](#)⁸¹

```
[C#]
public ref System.DateTime ILIMSQuery.FromDate { get; set; }
```

Parameters

Value

[out,retval]

2.3.7.5 ToDate Property

End date search range for object creation or modification date

```
[C/C++]  
HRESULT    get_ToDate( DATE* Value );  
HRESULT    put_ToDate( DATE Value );
```

```
[Visual Basic]  
Public Overloads Property ToDate As Date
```

VB code snippet See [ILIMSQuery::FromDate](#) ⁸³

```
[C#]  
public ref System.DateTime    ILIMSQuery.ToDate { get; set; }
```

Parameters

Value

[out,retval]

2.3.7.6 QueryDateType property

Type of date search, either a creation date, modification date, or both. See the enumerated type for LIMSQueryDateType.

```
[C/C++]  
HRESULT    get_QueryDateType( LIMSQueryDateType* Value );  
HRESULT    put_QueryDateType( LIMSQueryDateType Value );
```

```
[Visual Basic]  
Public Overloads Property QueryDateType As LIMSQueryDateType
```

VB code snippet See [ILIMSQuery::FromDate](#)⁸³

```
[C#]  
public LIMSQueryDateType ILIMSQuery.QueryDateType { get; set; }
```

Parameters

Value

[out,retval]

Note :

The qdtAllDateQuery value for the QueryDateType means that the date filters FromDate and ToDate are inactive. e.g the query will return all requested objects found.

2.3.7.7 ExecuteQuery Method

Execute the query

```
[C/C++]
HRESULT ExecuteQuery(
    ILIMSQueryResult** Result,
    ILIMSOperationResult** Value
);
```

```
[Visual Basic]
object.ExecuteQuery(
    ByRef Result As Object,
) As Object
```

VB code snippet See [ILIMSQuery::ObjectType](#)⁸¹

```
[C#]
void ILIMSQuery.ExecuteQuery(
    object Result,
);
```

Parameters

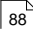
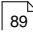
Result

[out]

Value

[out,retval]

2.3.8 ILIMSQueryResult

General methods and properties.		
property get	Count 	The number of objects returned by a query.
method	GetResultData 	method to get the individual results from a query

2.3.8.1 Count property

The number of objects returned by a query.

[C/C++]

```
HRESULT get_Count( LONG* Value );
```

[Visual Basic]

```
Public ReadOnly Property Count As Long
```

VB code snippet :

```
Public Function GetQueryResult()  
    Dim QIndx, Indx As Integer  
    Dim OpResult As LIMSOperationResult  
    Dim Qresult As LIMSTQueryResult  
    Dim Qdata As LIMSTQueryResultData  
    Dim LQuery As LIMSTQuery  
  
    If gConnection.LoggedIn = True Then  
        gConnection.Query.ObjectType = "Experiment" 'set the filters  
                                                SEE remarks below  
  
        gConnection.Query.Name = ""  
        gConnection.Query.Owner = ""  
        Set OpResult = gConnection.Query.ExecuteQuery(Qresult)  
        If OpResult.Successful = True Then  
            QIndx = Qresult.Count 'Get the number of items returned  
            If QIndx > 0 Then  
                Main.ExpList.Clear ' Delete the list in main frame  
                For Indx = 0 To QIndx - 1  
                    'Read Individual items and Add them to the list  
                    Set Qdata = Qresult.GetResultData(Indx)  
                    Main.ExpList.AddItem (Qdata.Name), Indx  
                Next Indx  
            End If  
        Else  
            MsgBox "Exp Query unsuccessful .."  
        End If  
    End If  
End Function
```

[C#]

```
public ref int ILIMSTQueryResult.Count { get; }
```

Parameters

Value

[out,retval]

2.3.8.2 GetResultData method

Method to get the individual results from a query.

```
[C/C++]
HRESULT GetResultData(
    LONG Index,
    ILIMSQueryResultData** Value
);
```

```
[Visual Basic]
object.GetResultData(
    ByVal Index As Long,
) As Object
```

VB code snippet : See [ILIMSQueryResult::Count](#) 

```
[C#]
void ILIMSQueryResult.GetResultData(
    int Index,
);
```

Parameters

Index

[in]

Value

[out,retval]

2.3.9 ILIMSQueryResultData

General methods and properties.		
property get	Name ⁹¹	Name of the object
property get	ObjectType ⁹¹	Type of the object
property get	CreationDate ⁹²	Creation date of the object
property get	ModificationDate ⁹²	Modification date of the object
property get	Path ⁹³	Folder path to the object on the database server

2.3.9.1 Name Property

Name of the object.

```
[C/C++]  
HRESULT    get_Name( BSTR* Value );
```

```
[Visual Basic]  
Public ReadOnly Property Name As String
```

```
[C#]  
public ref string    ILIMSQueryResultData.Name { get; }
```

Parameters

Value

[out,retval]

2.3.9.2 ObjectType Property

Type of the object.

```
[C/C++]  
HRESULT    get_ObjectType( BSTR* Value );
```

```
[Visual Basic]  
Public ReadOnly Property ObjectType As String
```

```
[C#]  
public ref string    ILIMSQueryResultData.ObjectType { get; }
```

Parameters

Value

[out,retval]

2.3.9.3 CreationDate Property

Creation date of the object.

```
[C/C++]  
HRESULT    get_CreationDate( DATE* Value );
```

```
[Visual Basic]  
Public ReadOnly Property CreationDate As Date
```

```
[C#]  
public ref System.DateTime    ILIMSQueryResultData.CreationDate { get; }
```

Parameters

Value

[out,retval]

2.3.9.4 ModificationDate Property

Modification date of the object.

```
[C/C++]  
HRESULT    get_ModificationDate( DATE* Value );
```

```
[Visual Basic]  
Public ReadOnly Property ModificationDate As Date
```

```
[C#]  
public ref System.DateTime    ILIMSQueryResultData.ModificationDate { get; }
```

Parameters

Value

[out,retval]

2.3.9.5 Path Property

Folder path to the object on the database server.

[C/C++]

```
HRESULT get_Path( BSTR* Value );
```

[Visual Basic]

```
Public ReadOnly Property Path As String
```

VB code snippet :

```
Public Function GetQueryResult()  
    Dim QIndx, Indx As Integer  
    Dim OpResult As LIMSOperationResult  
    Dim Qresult As LIMSQueryResult  
    Dim Qdata As LIMSQueryResultData  
    Dim LQuery As LIMSQuery  
    Dim ExpPath As String  
  
    If gConnection.LoggedIn = True Then  
        gConnection.Query.ObjectType = "Experiment" 'set the filters  
                                                    SEE ALSO ObjectType81  
  
        gConnection.Query.Name = ""  
        gConnection.Query.Owner = ""  
        Set OpResult = gConnection.Query.ExecuteQuery(Qresult)  
        If OpResult.Successful = True Then  
            QIndx = Qresult.Count 'Get the number of items returned  
            If QIndx > 0 Then  
                Main.ExpList.Clear ' Delete the list in main frame  
                For Indx = 0 To QIndx - 1  
                    'Read Individual items and Add them to the list  
                    Set Qdata = Qresult.GetResultData(Indx)  
                    ExpPath = Qdata.Path ' Get the path to experiment  
                    Main.ExpList.AddItem (Qdata.Name), Indx  
                Next Indx  
            End If  
        Else  
            MsgBox "Exp Query unsuccessful .."  
        End If  
    End If  
End Function
```

[C#]

```
public ref string ILIMSQueryResultData.Path { get; }
```

Parameters

Value

[out,retval]

2.4 Enumerated Types

This section contains information about the following enumerated types used with LIMSCClientLib.

- LIMSTypeQueryDateType⁹⁴

2.4.1 LIMSTypeQueryDateType

[C/C++]

```
typedef enum {  
    qdtModificationDateQuery = 1,  
    qdtCreationDateQuery     = 2,  
    qdtAllDateQuery          = 0  
} LIMSTypeQueryDateType;
```

[Visual Basic]

```
Enum LIMSTypeQueryDateType  
    qdtModificationDateQuery = 1  
    qdtCreationDateQuery     = 2  
    qdtAllDateQuery          = 0  
End Enum
```

[C#]

```
enum LIMSTypeQueryDateType {  
    qdtModificationDateQuery = 1,  
    qdtCreationDateQuery     = 2,  
    qdtAllDateQuery          = 0  
}
```

Constants

qdtModificationDateQuery

qdtCreationDateQuery

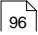
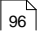
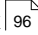
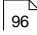
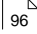
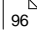
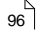
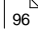
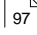
qdtAllDateQuery

Note :

The qdtAllDateQuery value for the QueryDateType means that the date filters FromDate and ToDate are inactive. E.g the query will return all requested objects found.

2.5 CoClasses

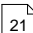
This section contains reference information for the COM CoClasses provided by LIMSCientLib. The following CoClasses are available with LIMSCientLib.

- LIMSConnection  96
- LIMSExperimentInfo  96
- LIMSInstrument  96
- LIMSOperationResult  96
- LIMSSampleDefinition  96
- LIMSSampleInfo  96
- LIMSQuery  96
- LIMSQueryResult  96
- LIMSQueryResultData  97

2.5.1 LIMSCONNECTION

Implements the ILIMSCONNECTION interface

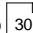
Implemented interfaces

- ILIMSCONNECTION 

2.5.2 LIMSEXPERIMENTINFO

Implements the ILIMSEXPERIMENTINFO interface

Implemented interfaces

- LIMSEXPERIMENTINFO 

2.5.3 LIMSinSTRUMENT

Implements the LIMSinSTRUMENT interface

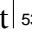
Implemented interfaces

- LIMSinSTRUMENT 

2.5.4 LIMSoPERATIONRESULT

Implements the LIMSoPERATIONRESULT interface

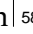
Implemented interfaces

- LIMSoPERATIONRESULT 

2.5.5 LIMSSAMPLEDEFINITION

Implements the LIMSSAMPLEDEFINITION interface

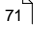
Implemented interfaces

- LIMSSAMPLEDEFINITION 

2.5.6 LIMSSAMPLEINFO

Implements the LIMSSAMPLEINFO interface

Implemented interfaces

- LIMSSAMPLEINFO 

2.5.7 LIMsQUERY

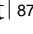
Implemented interfaces

- LIMsQUERY 

2.5.8 LIMsQUERYRESULT

Implements LIMsQUERYRESULT interface

Implemented interfaces

- LIMsQUERYRESULT 

2.5.9 LIMSQueryResultData

Implements ILIMSQueryResultData

Implemented interfaces

- ILIMSQueryResultData

3 XML output

This chapter describes the complete format of the XML output, as returned by the LIMS interface when an experiment summary is requested.

3.1 Additional details, requirements

specifying the results returned by LIMS

An Abs Quant result will include the following for each sample

- Cp call
- Calculated concentration
- Standard concentration
- Interpolated/Extrapolated flag for concentration
- Certain/uncertain flag on the Cp call
- High/Low/Normal flag for the Cp call

In the Tm module, the results will include a flag to indicate manual editing for each sample.

➤ For each sample and for each of the six peaks that may be called in the Tm algorithm, the Tm calling results will include the following:

- Tm
- Area
- Peak
- Width
- Shoulders

➤ Rel Quant results will include the following for result set

- Result set name (this column is empty for calibrators)
- Sample Type, position, and sample name for samples in the results set. Values for sample type are Target Calibrator, Reference Calibrator, Target Unknown, and Reference Unknown.
- Cp call for each sample in the result set and the median Cp for the set.
- Concentration Ratio and error (Concentration ratio of the Calibrators or of the Target and Reference Unknowns in this result set)
- Normalized Ratio and error (The normalized ratio of the Unknowns with the Calibrators for this result set)
- Multiplication/Correction Factor (The multiplication and correction factors for the result set, expressed as a fraction.)

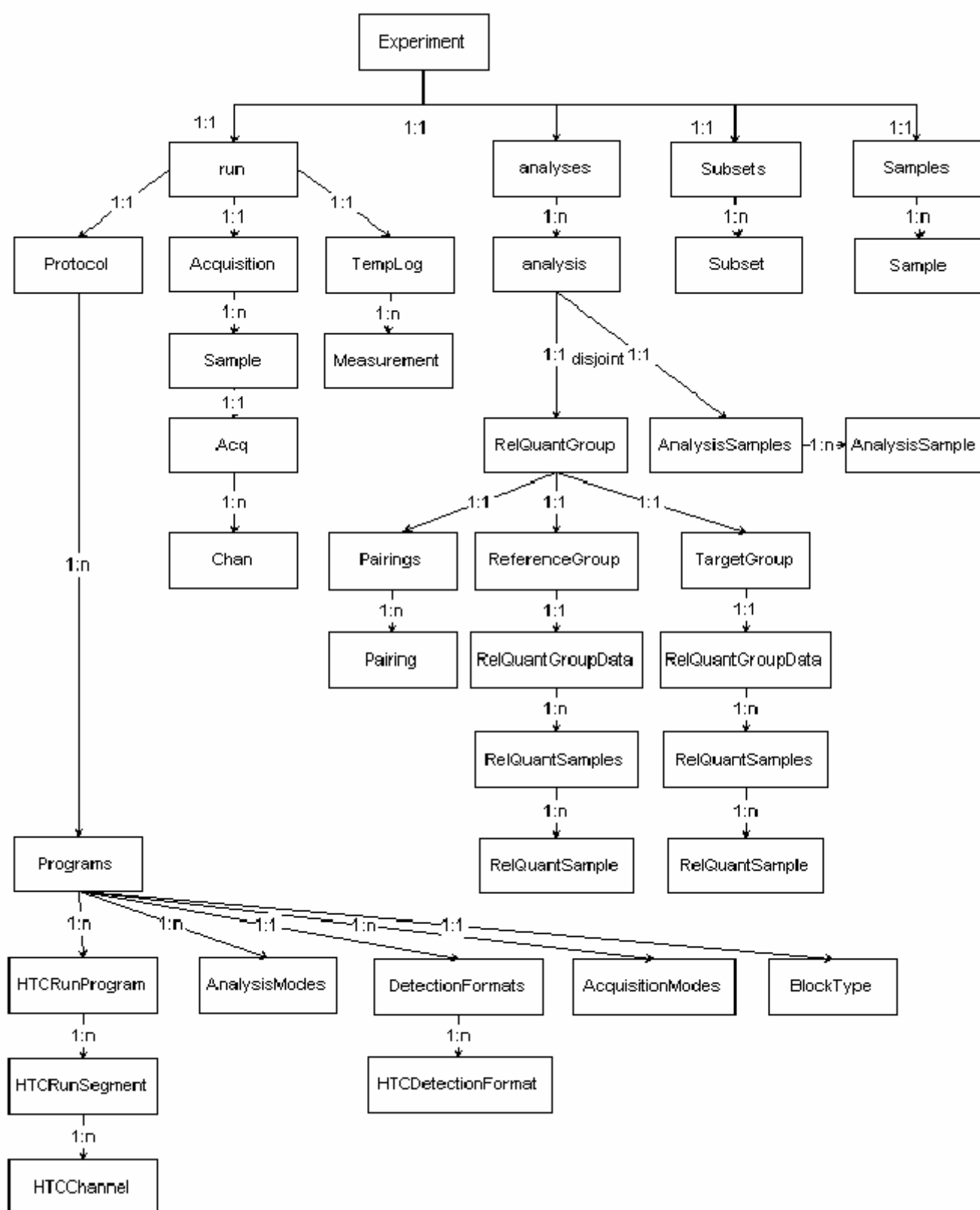
➤ Genotyping results from LIMS will include the following for each sample:

- Group Name
- Score
- Resolution
- A flag that indicates manual editing

➤ LIMS shall return the following for each sample in a Fit Points analysis:

- Cp call
- Calculated concentration
- Standard concentration
- Interpolated/Extrapolated flag for call and concentration

3.2 Output file schema



3.3 XML-File Elements

3.3.1 Experiment

Experiment is the root element of the XML-file.

- Child Elements:

run
analyses
Subsets
Samples

- Attributes:

Name	String	
Created	Date	yyyy-mm-ddThh:mm:ss.msc
createdByName	String	
LastModified	Date	yyyy-mm-ddThh:mm:ss.msc
LastModifiedByName	String	
SWVersion	String	LCS480 1.0.0.95
RevsComplete	Integer	>0

3.3.1.1 Run

- Child Elements:

Protocol
Acquisition
TempLog

- Attributes:

Name	String	
Created	Date	yyyy-mm-ddThh:mm:ss.msc
CreatedByName	String	
LastModified	Date	yyyy-mm-ddThh:mm:ss.msc
State	String	rscompleted
StartTime	Date	yyyy-mm-ddThh:mm:ss.msc
EndTime	Date	yyyy-mm-ddThh:mm:ss.msc
InstrumentID	Integer	
InstrumentVersion	String	HTC_VER_10
InstrumentName	String	Pilot 517

InstrCalibrationDate	Date	yyyy-mm-dd hh:mm:ss:ms
Technician	String	
Notes	Text	

3.3.1.1.1 Protocol

- Child Elements:

Programs

- Attributes:

Name	Type	Format
class	String	
version	Integer	
created	Date	yyyy-mm-ddThh:mm:ss.msc
last modified	Date	yyyy-mm-ddThh:mm:ss.msc

- Child Elements:

Emlist	list of HTCRunProgram nodes
AnalysisModes	list of Strings
DetectionFormats	Node with list of HTCDetectionFormat
AcquisitionModes	List of Strings
BlockType	

- Attributes:

class	String	
version	Integer	
ChannelCount	Integer	
InstrumentSubclass	String	
SeekTemp	Integer	
MaxPositionsToSeek	Integer	
SampleVolume	Integer	

re.

- Child Elements:

HTCRunSegment

- Attributes:

version	Integer	
Name	String	
Cycles	Integer	Pre-incubation, Amplification
AnalysisMode	String	Quantification

.

- **Attributes:**

Hold	Integer	
Slope	Float	
StepDelay	Integer	
StepSize	Integer	
Target	Integer	
Target2	Integer	
AcquisitionMode	Integer	
AcqPerDegree	Integer	

- **Child Elements:**

HTCDetectionFormat	
--------------------	--

- **Attributes:**

class	String	
version	Integer	
DefFormatNdx	Integer	

- **Attributes:**

class	String
version	Integer
Id	String
RowCount	Integer
ColCount	Integer
OvershootDnDelay	Integer
OvershootUpDelay	Integer

RampRateMaxDn	Float
RampRateMaxUp	Float
ReactionVolMin	Integer
ReactionVolMax	Integer
ReactionVolDefault	Integer
HorVertCrosstalkCoefficient	Float
DiagCrosstalkCoefficient	Float

3.3.1.1.2 Acquisition

3.3.1.1.3 Acquisition

- Child Elements:

Acquisition with child element Sample
Sample with child element Acq
Acq with child element Chan

- Attributes of Chan:

Fluor	float	
Temp	float	
Time	integer	>0

3.3.1.1.4 TempLog

- Child Elements:

Measurement

- Attributes:

Temp	float	
Time	integer	>0

3.3.1.2 Analyses

3.3.1.2.1 Analysis

- Elements:

Pairings	Mapping of ReferenceGroup to TargetGroup
ReferenceGroup	
TargetGroup	

Elements:

RelQuantSamples

.

- Attributes:

name	String	
Position	String	A1, D5
RQSampleIncluded	Boolean	0,1
RQSampleType	String	
RQSampleCrossingPoint	Float	
RQSampleCall	String	pdCPositive
RQSampleConcentration	Float	
Incomplete	integer	
name	String	
Position	String	A1, D5
RQSampleIncluded	Boolean	0,1
RQSampleType	String	

- Elements:

Pairing

- Attributes:

Name	String	
------	--------	--

RQResultSetNormalizedRatio	Float	
RQResultSetConcentrationRatio	Float	
RQResultSetCalibratorConcentrationRatio	Float	
RQResultSetState	String	rsPositive
RQResultSetNormalizedRatioError	Float	
RQResultSetConcentrationRatioError	Float	
RQResultSetCalibratorConcentrationRatioError	Float	
RQResultSetMultFactor	Integer	
RQResultSetCorrFactor	Integer	
CalRef	Integer	
CalTarget	Integer	
Ref	Integer	
Target	Integer	
Target Median Cp	Float	
Reference Median Cp	Float	

Attributes:

name	String	
Position	String	A1, D5
IsIncluded	Boolean	0,1
CrossingPoint	Float	
CpUncertain	Boolean	0,1
CpState	String	
CalcConc	Float	
StandardConc	Float	
CalcConcUnc	Boolean	0,1

.

- Attributes:

GroupName
Res
Score
ManualGroup

- Attributes:

name	String	
Position	String	A1, D5
IsIncluded	Boolean	0,1
Call	Integer	
TmCount	Integer	
ManualTms	Integer	

- Lists:

Shoulder	Array of float	TM of shoulder if any occur
Tms	Array of float	number of peaks
Amounts	Array of float	value of the area of the peak
Heights	Array of float	height of the peak
Widths	Array of float	width of the peak

Attributes:

name	String	
Position	String	A1, D5
IsIncluded	Boolean	0,1

.

3.3.1.3 Subsets

List of subsets

3.3.1.4 Samples

List of samples

Index

- C -

collection 58, 59
connection object 21
Control unit 16

- D -

database 19

- E -

Environment 19
Exor 10, 16
Exor DB 12
ExperimentInfo 24

- H -

heartbeat 12

- I -

Instrument 16
Interface 21
IXO file 35

- L -

Language settings 19

- M -

macro 10, 42
Micro Well Plate 10
MWP 10

- O -

Operating System 19

- R -

Restrictions 19

- S -

samples
 Sample collection 58
synchronous 12

- T -

time stamp 56

- U -

user 18
User parameter 28

This page is left for your personal notes. □

