

# Chapter 4 Data Transmission (LIS connection)

## 4.1 Function of LIS connection:

Transmit the urinary sediment result to the LIS side, and then print out the report sheet with the LIS software.

## 4.2 LIS connection operation:

1. The default LIS route is D:\LISText, and modification can also be made in setup>system screen, as shown in Fig. 4-1.

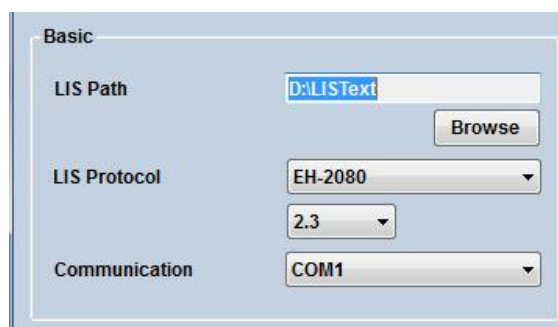


Fig. 4-1 Screen for setting up system module

2. After completing specimen result review, click “LIS key” on the software Review screen or report screen, the exported text file can be seen in the LIS export route (D:\LISText) set in step 1, and for each specimen, a text file will be generated.



Fig. 4-2 Review and report screens

**Note: In connecting the LIS, it is only necessary to tell LIS engineers this route, and they will set up connection according to the text message.**