


USER INTERFACE RS232-C

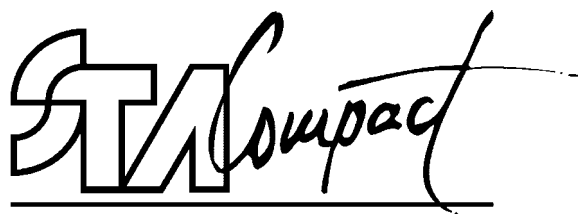

Cat. No.: 28062

Revision Table

Manual Version	Date	List of Modifications
V 1.0	November 95	Original Version
V 1.1.	March 97	Modified Pages: 11, 14, 24 and Table of Contents Note: Modifications made to be consistent with software version 104.
V 1.1	March 02	Addition of an automatically incremented sequence number (in the frame Order) transmitted to the host

1.	Physical Characteristics	1
2.	Standards	1
3.	Terminology	1
3.1.	Global Definitions	1
3.2.	Definition of Characters	2
4.	Global Options Menu	3
5.	Exchange Protocol	4
5.1.	Definition of Line Modes on the STA Compact [®]	4
5.2.	Definition of Phases of Communication	4
5.3.	Host Computer → STA Compact [®] Exchanges	4
5.4.	STA Compact [®] → Host Computer Exchanges	5
6.	Record Definitions	6
6.1.	Message Header Record	6
6.2.	Patient Information Record	7
6.3.	Test Order Record	8
6.4.	Result Record	10
6.5.	Request Information Record	12
6.6.	Message Terminator Record	12
6.7.	Manufacturer Information Record	13
7.	Typical Contents of Some Messages	15
7.1.	Result Transmission (STA Compact [®] → Host)	15
7.2.	Request for a Working List	15
7.3.	Working List Return	15
8.	Flowcharts (ASTM)	16
8.1.	Sequencing of Reception/Transmission	16
8.2.	Acquisition of a Record from the Host Computer	17
8.3.	Establishment Phase (STA Compact [®] → Host Computer)	18

8.4	Transfer Phase (STA Compact[®] → Host Computer)	19
8.5.	Termination Phase (STA Compact[®] → Host Computer)	20
9.	Examples of Traces (ASTM)	21
9.1.	Transmission of Patient File (STA Compact[®] → Host Computer)	21
9.2.	Request for a Working List	23
9.3	Working List Return	24
9.4.	Transmission of a Quality Control Result (STA Compact[®] → Host Computer)	25
9.5.	On-Line Test	25

1. Physical Characteristics

- Interface: Series type, RS232C standard
- Connections: DB-9 male plug on STA Compact[®]
pin 2: reception input
pin 3: transmission output
pin 5: ground
- Exchange protocol: no hardware control line is connected;
- Baud Rate^{*}: 300, 600, 1200, 2400, 4800 or 9600 baud;
- Parity^{*}: even, odd or none;
- Data Bit^{*}: 7 or 8;
- Stop Bit^{*}: 1 or 2.

^{*}: parameters to be programmed in the **Global Options** menu in the STA Compact[®] software (see below, chapter 4.4.4. of this STA Compact[®] User Interface RS232-C).

2. Standards

The ASTM connection protocol on STA Compact[®] is based on the two following standards:

- E 1381-94 : Specification for Low-Level Protocol to transfer messages between Clinical Laboratory Instruments and Computer Systems,
- E 1394-91 : Standard Specification for transferring information between Clinical Instruments and Computer Systems.

3. Terminology

3.1. Global Definitions

Sender:	The device that has a message to send and that initiates transmission process.
Upload:	Data transmitted from the instrument (STA Compact [®]) to the host computer.
Receiver:	The device that responds to the sender.
Download:	Data transmitted from the host computer to the instrument (STA Compact [®]).
Frame:	All framing required for transmission of data. This framing will include: <STX> [Frame Number] [DATA] <ETB> or <ETX> [C1] [C2] <CR> <LF>.
Frame Number: (Fn)	Frame number permits the receiver to distinguish between new and retransmitted frames. ASCII digit ranging from 0 to 7. The frame number begin at 1 with the first frame of transfer phase. It is incremented by one for every new frame transmitted. After 7, the frame number rolls over to 0, and continue in this fashion.

Message: A group of records that begins with a Header record, and ends with a Terminator record.

Record: This is the DATA as described within the frame. If the DATA is longer than 240 characters, then it must be split into two parts and sent in two frames. The intermediate frame will use the <ETB> character, and the ending frame will use the <ETX> character.

For our application:

- in download mode, a frame may include several records, indeed a message;
- in upload mode, one frame is used for sending each record.

List of records selected:

- H (Header);
- P (Patient);
- O (Order);
- R (Result);
- Q (Query);
- L (Terminator);
- M (Manufacturer).

(See description, Chapter 6. of this STA Compact[®] User Interface RS232-C).

3.2. Definition of Characters

<ACK> (06h)	:	acknowledge without error;
<NACK> (15h)	:	acknowledge with error;
<ETX> (03h)	:	end of text transmission;
<STX> (02h)	:	start of text transmission;
[C1]	:	most significant character of checksum;
[C2]	:	least significant character of checksum;
<ENQ> (05h)	:	request for line (enquiry);
<ETB> (17h)	:	end of transmission block;
<EOT> (04h)	:	end of transmission;
<CR> (0Dh)	:	carriage return;
<LF> (0Ah)	:	line feed.

Calculation of Frame Checksum = the checksum is computed by adding the binary values of the characters (except for <STX>, [C1], [C2], <CR> and <LF>), keeping the least significant eight bits of the result. The checksum is encoded as two characters. (For example a checksum of 122 can be represented as 7A in hexadecimal. The checksum is transmitted as the ASCII character 7 followed by the character A.).

4. Global Options Menu

Access to the Global Options menu:

From the main screen:

- press the **ESC** key to bring up the **Main Menu**;
- select the **Setup** menu, then select the **Global Options** sub-menu from this menu.

This menu in the STA Compact® program is used, amongst other things, for programming the communication parameters between an STA Compact® and the host computer. Access is restricted by a code system.

Fig. 1 - GLOBAL OPTIONS screen, page 2 of 2 (ASTM Protocol)

Esc Quit
PgUp Page 1

GLOBAL OPTIONS
Page 2/2

COMMUNICATIONS		FILE ACQUISITION	
		Item	Format
Protocol ASTM		1	16 characters
Station Number (0..99)	99	2	16 -
Baud Rate (bauds)	9600	3	12 -
Parity	none	4	6 -
Number of Data Bits	8	5	4 -
Number of Stop Bits	1	ID Type	
Number of On Error Retries	3	AlphaNum	
Verify Patient Data	no	MISCELLANEOUS	
Send Sequence Number	no	Time Format	
		Date Format	

Description of Parameters

Parameters for which there is a multiple choice are shown in white (to make the different selections, use the **DELETE** key), and those with free entry in blue.

- *Station Number*: number by which the central site identifies the STA Compact®; Maximum: 99;
- *Baud Rate*: chosen transmission speed: 300, 600, 1200, 2400, 4800 or 9600 baud;
- *Parity*: type of parity chosen: none, even or odd;
- *Number of Data Bits*: 7 or 8;
- *Number of Stop Bits*: 1 or 2;
- *Number of On Error Retries*: number of times the STA Compact® returns the frame when it receives a character different from <ACK> or <EOT>. In accordance with the Standard E1381-91, this number has to be set to 6.

- *Verify Patient data:*

possible choices:

- **YES**, when the working list is remotely loaded, the STA Compact[®] records 4 patient information fields (a maximum of 38 characters for processing and control) and the operator may check them. These 38 characters are defined on the STA Compact[®] as follows:

Information 1	16 characters maximum
Information 2	12 characters maximum
Information 3	6 characters maximum
Information 4	4 characters maximum
Total	38 characters maximum

- **NO**, when the working list is remotely loaded, the STA Compact[®] records the 4 patient information fields but the operator cannot check them.

For more the details, see *Chapter 6.6.2. of this STA Compact[®] User Interface RS232-C*.

- *Send Sequence Number*: automatically incremented number generation (between 1 through 9999, number reset at 9999) located in the **Order** frame, field 4 that is the **Instrument Specimen ID** field (not activated by default).

5. Exchange Protocol

5.1. Definition of Line Modes on the STA Compact[®]

- Idle = no communication in progress
- Tx = transmission by STA Compact[®]
- Rx = reception by STA Compact[®]

5.2. Definition of Phases of Communication

ASTM communication is in three stages:

- Establishment phase: line picked up by the host computer or the STA Compact[®];
- Transfer phase: transmission of the different records making up the message (*see definition of records, chapter 6. of the STA Compact[®] User Interface RS232-C*);
- Termination phase: the line is released.

5.3. Host Computer → STA Compact[®] Exchanges

STA Compact[®] responses to characters transmitted by the host computer:

1. Reception of the <ENQ> character transmitted by the host computer:
If the STA Compact[®] is in the **Idle** line mode, it transmits an <ACK> character (06h) and then changes to reception line mode **Rx**.
If after 30 seconds, no character is received, then STA Compact[®] returns to the **Idle** line mode.
2. Reception of the <EOT> character transmitted by the host computer:
 - If the STA Compact[®] is in the **Rx** reception line mode and a header record has been received, the STA Compact[®] then changes to **Idle** line mode (releasing the line); the STA Compact[®] then analyses the records received;
 - If the STA Compact[®] is in **Tx** transmission line mode, the <EOT> character is regarded as an acknowledgement (*see chapter 5.4. b. of the STA Compact[®] User Interface RS232-C*).
3. Reception of other characters transmitted by the host computer:
The STA Compact[®] proceeds to enter the characters received. After entering a complete frame, the STA Compact[®] verifies the checksum. If this is satisfactory, the STA Compact[®] acknowledges by returning an <ACK> character (06h). If not, a <NACK> character (15h) is returned.
Following the upload of the characters <ACK> or <NACK>, if after 30 seconds, no character is received, then STA Compact[®] returns to the **Idle** line mode.

5.4. STA Compact[®] → Host Computer Exchanges

a) Establishment Phase

If the line mode is **Idle**, the STA Compact[®] sends an <ENQ> character (enquiry).

The following are the possible responses from the host computer:

1. Transmission of an <ACK> Character by the Host Computer:
The STA Compact[®] changes to **Tx** transmission line mode and initiates the record transfer phase if it is in normal mode or the termination phase (transmission of an <EOT> character) if it is in line test mode.
2. Transmission of an <ENQ> Character by the Host Computer:
The STA Compact[®] waits 5 seconds and then returns an <ENQ> character.
3. Transmission of Other Characters by the Host Computer:
The STA Compact[®] waits 10 seconds and then returns the <ENQ> character. When the value of the "Number of On Error Retries" defined in the **Global Options** menu (see chapter 4. of the *STA Compact[®] User Interface RS232-C*) is reached, the STA Compact[®] changes to termination phase (transmission of an <EOT> character) and displays an error message.
4. No Response from the Host Computer:
The STA Compact[®] waits 15 seconds, then changes to termination phase (transmission of an <EOT> character) and displays an error message.

b) Transfer Phase (STA Compact[®] → Host Computer)

Transmission of the first record by the STA Compact[®].

The following are the possible responses of the host computer:

1. Transmission of the <ACK> or <EOT> Character by the Host Computer:
The STA Compact[®] sends the next record.



Caution! The STA Compact[®] regards the <EOT> character as an acknowledgement character and continues to transmit records.

2. Transmission of Other Characters by the Host Computer:
The STA Compact[®] waits 10 seconds before returning the record. If the number of on Error Retries is higher than that defined in the **Global Options** menu (see chapter 4. of the *STA Compact[®] User Interface RS232-C*), the STA Compact[®] changes to termination phase (transmission of an <EOT> character) and displays an error message.
3. No Response from the Host Computer:
The STA Compact[®] waits 15 seconds and then changes to termination phase (transmission of an <EOT> character) and displays an error message.

c) Termination Phase (STA Compact[®] → Host Computer)

The STA Compact[®] transmits an <EOT> character and does not wait for any return character. The STA Compact[®] changes to **Idle** line mode.

d) Example of Host Computer → STA Compact® Exchange

STA Compact®	Host Computer
<ENQ> ----->	
<STX> 1 [DATA1] <ETX> [C1] [C2] <CR> <LF>	<ACK> or <EOT>
<STX> 2 [DATA2] <ETX> [C1] [C2] <CR> <LF>	<ACK> or <EOT>
<STX> 3 [DATA3] <ETX> [C1] [C2] <CR> <LF>	<ACK> or <EOT>
<EOT> ----->	No response expected

6. Record Definitions**6.1. Message Header Record**

Purpose : To identify the sender.

STA Compact® implementation example:

H | \ ^ & | | | 99 ^ 2.00 | | | | | P | 1.00 | 19941213100400 <CR>

Description of a Message Header Record on STA Compact®:

PARAGRAPH'S NUMBER OF E 1394-91 STANDARD	CHARACTER	FIELD DESCRIPTION	SUPPORTED	
			UPLOAD (STA COMPACT® → HOST)	DOWNLOAD (HOST → STA COMPACT®)
7.1.1	H	Record Type ID = header record	Yes	Yes
7.1.2	 \ ^ &	Delimiter Definition = defined as field delimiter (7Ch) = defined as repeat field delimiter (5Ch) = defined as component delimiter (5Eh) = defined as Escape delimiter (26h)	Yes	Yes
7.1.3		Message Control ID	No	No
7.1.4		Access Password	No	No
7.1.5	 99 ^ 2.00	Sender Name or ID = Instrument ID ❶ (Numeric 1 to 99) = Version level of Stago instrument implementation	Yes Yes	Yes Yes
7.1.6		Sender Street Address	No	No
7.1.7		Reserved Field	No	No
7.1.8		Sender Telephone Number	No	No

PARAGRAPH'S NUMBER OF E 1394-91 STANDARD	CHARACTER	FIELD DESCRIPTION	SUPPORTED	
			UPLOAD (STA COMPACT [®] → HOST)	DOWNLOAD (HOST → STA COMPACT [®])
7.1.9		Characteristics of sender	No	No
7.1.10		Receiver ID	No	No
7.1.11		Comment or Special Instructions	No	No
7.1.12	 P	Processing ID P = Production Q = Quality Control T = Training D = Debugging	Yes Yes No No	No No No No
7.1.13	 1.00	Version No. = Version level of ASTM specifications (E 1394-91)	Yes	No
7.1.14	 19941213 100400	Date and Time of Message Format: yyyyymmddhhmmss yyyy = year mm = month dd = day hh = hour mm = minutes ss = seconds	Yes	No
	<CR>	carriage return		

❶ : For STA Compact[®], this parameter is defined within the **Global Options** menu (see Chapter 4. of this STA Compact[®] User Interface RS232-C).

6.2. Patient Information Record

Purpose : To give information on the Patient.

STA Compact[®] implementation example:

P | 1 | | | Info1^Info2^Info3^Info4 <CR>

Description of a Patient Information Record on STA Compact[®]:

PARAGRAPH'S NUMBER OF E 1394-91 STANDARD	CHARACTER	FIELD DESCRIPTION	SUPPORTED	
			UPLOAD (STA COMPACT [®] → HOST)	DOWNLOAD (HOST → STA COMPACT [®])
8.1.1	P	Record Type ID = Patient record (50h)	Yes	Yes
8.1.2	 1	Sequence Number	Yes	Yes
8.1.3		Practice Assigned Patient ID	No	No
8.1.4		Laboratory Assigned Patient ID	No	No

PARAGRAPH'S NUMBER OF E 1394-91 STANDARD	CHARACTER	FIELD DESCRIPTION	SUPPORTED	
			UPLOAD (STA COMPACT [®] → HOST)	DOWNLOAD (HOST → STA COMPACT [®])
8.1.5		Patient ID No.3	Yes	Yes
	Info1	= 16 Alphanumeric Characters Max.		
	^			
	Info2	= 12 Alphanumeric Characters Max.		
	^			
	Info3	= 6 Alphanumeric Characters Max.		
	^			
	Info4	= 4 Alphanumeric Characters Max.		
	<CR>	Carriage Return		

Note: Fields 8.1.6 to 8.1.35 as defined in the ASTM Specifications E 1394-91 are not supported. In downloading (Host → STA Compact[®]), if these fields exist, they are ignored. In uploading (STA Compact[®] → Host), we use the short record format described above.

6.3. Test Order Record

Purpose:

- In downloading (Host → STA Compact[®]), this record is used by the host computer to give orders to the instrument (STA Compact[®]): list of tests to be run;
- In uploading (STA Compact[®] → Host), this record is used by STA Compact[®] to give the sample identification.

STA Compact[®] implementation example for download (Host → STA Compact[®]):

O | 1 | 00120025 | | ^10 \ ^11 \ ^12 | R <CR>

STA Compact[®] implementation example for upload (STA Compact[®] → Host):

O | 1 | 00120025 | | | R <CR>

Description of a Test Order Record:

PARAGRAPH'S NUMBER OF E 1394-91 STANDARD	CHARACTER	FIELD DESCRIPTION	SUPPORTED	
			UPLOAD (STA COMPACT [®] → HOST)	DOWNLOAD (HOST → STA COMPACT [®])
9.4.1	O	Record Type ID = Order record	Yes	Yes
9.4.2	 1	Sequence Number	Yes	Yes
9.4.3	 00120025	Specimen ID = 16 Alphanumeric Characters Max.	Yes	Yes

PARAGRAPH'S NUMBER OF E 1394-91 STANDARD	CHARACTER	FIELD DESCRIPTION	SUPPORTED	
			UPLOAD (STA COMPACT® → HOST)	DOWNLOAD (HOST → STA COMPACT®)
9.4.4		Instrument Specimen ID <i>Send Sequence Number</i> (in the Global Options menu): ⇒ activated ⇒ not activated	Yes No	No No
9.4.5	 ^ ^ ^ 10 \ ^ ^ ^ 11 \ ^ ^ ^ 12	Universal Test ID Universal Test ID Part 1 Universal Test ID Name Universal Test ID Type Manufacturer's or Local Code = ID of the first requested test ❶ = Repeat field delimiter = ID of the second requested test ❶ = Repeat field delimiter = ID of the third requested test ❶	No	Yes
9.4.6	 R	Priority R = Routine S = Stat A = As soon as possible C = Callback P = Preoperative	Yes Yes No No No	Yes Yes No No No
	<CR>	Carriage Return		

❶ For STA Compact®, these identifications are the transmission test numbers defined for each test main units (see *STA Compact® Operator's Manual, Chapter 3.1.3.1.*). For each sample, STA Compact® can accept a maximum of **twelve tests**.

Note: Field 9.4.7 to 9.4.31 as defined in the ASTM Specifications E 1394-91 are not supported. In downloading (Host → STA Compact®), if these fields exist, they are ignored. In uploading (STA Compact® → Host), we use the short record format described above.

6.4. Result Record

Purpose: To transmit Results to the host computer.

STA Compact® implementation example:

R | 1 | ^10 | 100 | % | | | F | | | 19941213100400 <CR>

Description of a Result Record on STA Compact®:

PARAGRAPH'S NUMBER OF E 1394-91 STANDARD	CHARACTER	FIELD DESCRIPTION	SUPPORTED	
			UPLOAD (STA COMPACT® → HOST)	DOWNLOAD (HOST → STA COMPACT®)
10.1.1	R	Record Type ID = Result Record (52h)	Yes	No
10.1.2	 1	Sequence Number	Yes	No
10.1.3	 ^ ^ ^ 10	Universal Test ID Universal Test ID Part 1 Universal Test ID Name Universal Test ID Type Manufacturer's or Local Code = ID of the test ❶	Yes	No
10.1.4	 100	Data or Measurement Value = Test result	Yes	No
10.1.5	 %	Units = Unit of the test result ❷	Yes	No
10.1.6		Reference Ranges	No	No
10.1.7		Result Abnormal Flags	No	No
10.1.8		Nature of Abnormality Testing	No	No
10.1.9	 F	Result Status F = Final results C = Correction of previously transmitted results P = Preliminary results X = Results cannot be done, request will not be honored I = In instrument, results pending S = Partial results M = This result is a MIC level R = This result was previously transmitted	Yes No No No No No No No	No No No No No No No No

PARAGRAPH'S NUMBER OF E 1394-91 STANDARD	CHARACTER	FIELD DESCRIPTION	SUPPORTED	
			UPLOAD (STA COMPACT [®] → HOST)	DOWNLOAD (HOST → STA COMPACT [®])
		N = This result record contains necessary information to run a new order Q = This result is a response to an outstanding query V = Operator verified/approved result	No	No
10.1.10		Date of Change in Instrument Normative Values or Units	No	No
10.1.11		Operator Identification	No	No
10.1.12		Date/Time Test Started	No	No
10.1.13	 19941213 100400	Date/Time Test Completed Format: yyyymmddhhmmss yyyy = year mm = month dd = day hh = hour mm = minutes ss = seconds	Yes ③	No
	<CR>	Carriage return		



Caution! One record is used for each result. If for a requested test, result can be expressed in 3 different units, then there will be 3 records.

- ① For STA Compact[®], this identification is the transmission test number defined for each related test unit (see *STA Compact[®] Operator's Manual, Chapter 3.1.3.1.*)
- ② List of units available on STA Compact[®]:

Units	Result Format
sec., Reference Time	xxx.x
%, mg/dl	xxxxx
OD, OD/min	x.xxx
Ratio, INR, g/l, ng/ml, U/ml, IU/ml, µg/ml	xx.xx
UA1, UA2, UA3	0 to 3 decimal

- ③ Is only used while uploading quality control results.



Caution! On STA Compact[®], a Result Record is **always followed** by a Manufacturer Information Record (see *below, Chapter 6.7.*).

Note: Field 10.1.14 as defined in the ASTM Specifications E 1394-91 is not supported. In downloading (Host → STA Compact[®]), if this field exist, it is ignored. In uploading (STA Compact[®] → Host), we use the short record format described above.

6.5. Request Information Record

Purpose: For instrument (STA Compact[®]), to request the working list (tests to be run).

STA Compact[®] implementation example:

Q | 1 | ^ 00120025 <CR>

Description of a Request Information Record:

PARAGRAPH'S NUMBER OF E 1394-91 STANDARD	CHARACTER	FIELD DESCRIPTION	SUPPORTED	
			UPLOAD (STA COMPACT [®] → HOST)	DOWNLOAD (HOST → STA COMPACT [®])
12.1.1	Q	Record Type ID = Request record	Yes	No
12.1.2	 1	Sequence Number	Yes	No
12.1.3	 ^ 00120025	Starting Range ID Number Computer System Patient ID Speciment ID = 16 Alphanumeric Characters Max.	Yes	No
	<CR>	Carriage return		

Note: Field 12.1.4 to 12.1.13 as defined in the ASTM Specifications E 1394-91 are not supported. In downloading (Host → STA Compact[®]), if these fields exist, they are ignored. In uplaoding (STA Compact[®] → Host), we use the short record format described above.

6.6. Message Terminator Record

Purpose: Last record of message.

STA Compact[®] implementation example:

L | 1 | N <CR>

Description of a Terminator Record:

PARAGRAPH'S NUMBER OF E 1394-91 STANDARD	CHARACTER	FIELD DESCRIPTION	SUPPORTED	
			UPLOAD (STA COMPACT [®] → HOST)	DOWNLOAD (HOST → STA COMPACT [®])
13.1.1	L	Record Type ID = Terminator record	Yes	Yes
13.1.2	 1	Sequence Number	Yes	Yes

PARAGRAPH'S NUMBER OF E 1394-91 STANDARD	CHARACTER	FIELD DESCRIPTION	SUPPORTED	
			UPLOAD (STA COMPACT [®] → HOST)	DOWNLOAD (HOST → STA COMPACT [®])
13.1.3	 N	Terminator Code N = Normal termination T = Sender aborted R = Receiver requested abort E = Unknow system error Q = Error in last request for information I = No information available from last query F = Last request for infor- mation processed	Yes No No No No No No	No No No No No No No
	<CR>	Carriage return		

6.7. Manufacturer Information Record

Purpose: To transmit result error code and result alarm code.

STA Compact[®] implementation example:

M | 1 | A | @ <CR>

Description of a Manufacturer Information Record:

PARAGRAPH'S NUMBER OF E 1394-91 STANDARD	CHARACTER	FIELD DESCRIPTION	SUPPORTED	
			UPLOAD (STA COMPACT [®] → HOST)	DOWNLOAD (HOST → STA COMPACT [®])
15.1.1	M	Record Type ID = Manufacturer Information Record	Yes	No
15.1.2	 1	Sequence Number	Yes	No
15.1.3	 A	Error Code A = Validated (see error code list below)	Yes	No
15.1.4	 @	Alarm Code @ = No alarm (see alarm code list below)	Yes	No
	<CR>	Carriage return		



Caution! On STA Compact[®], a Manufacturer Information Record **always follows** a Result Record to indicate the Error and the Alarm Codes of the result.

Examples:

R | 1 | ^10 | 100 | % | | | F | | | 19941213100400 <CR>

M | 1 | A | @ <CR>

R | 2 | ^11 | 13.2 | sec. | | | F | | | 19941213100400 <CR>

M | 2 | A | @ <CR>

R | 3 | ^12 | 1.23 | INR | | | F | | | 19941213100400 <CR>

M | 3 | A | @ <CR>

.....

Definition of Error Codes:

Characters	ASCII Code	Error Code Definition
'1'	31h	"to be validated"
'2'	32h	"tech error"
'3'	33h	"> Mmax."
'4'	34h	"< Mmin."
'5'	35h	"diff > tol"
'6'	36h	"QNS", no plasma
'8'	38h	"linearity"
'A'	41h	"Validated"

Definition of Alarm Codes:

Characters	ASCII Code	Alarm Description
'@'	40h	No alarm
'A'	41h	Result : Confirmed with T>max.
'B'	42h	Not used
'C'	43h	Quality Control : out of range or not done
'D'	44h	Quality Control : overridden
'E'	45h	Needle #3 : no level detection
'F'	46h	Needle #2 : no level detection
'G'	47h	Needle #1 : no level detection
'H'	48h	Result : value in primary units skewed
'I'	49h	Result : dilution change
'J'	4Ah	Result : rerun test
'K'	4Bh	Reagent drawer : temperature out of limit
'L'	4Ch	Syringe : maintenance date overdue
'M'	4Dh	Not used
'N'	4Eh	Not used

Note: If there are several alarm codes for a result, then the most priority one is transmitted.
Alarm codes are classified from '@' to 'N' with highest priority for code '@'.

7. Typical Contents of Some Messages

7.1. Result Transmission (STA Compact[®] → Host)

Each message refers to only one sample:

- Header record;
- Patient record;
- Test Order record;
- Result record;
- Manufacturer Information record;
- Message Terminator record.

7.2. Request for a Working List

Only one sample is dealt with per message:

- Header record;
- Request Information record;
- Message Terminator record.

7.3. Working List Return


Only one sample is dealt with per message:

- Header record;
- Patient record;
- Test Order record;
- Message Terminator record.

8. Flowcharts (ASTM)


8.1. Sequencing of Reception/Transmission

Fig. 2 - Sequencing of Reception/Transmission

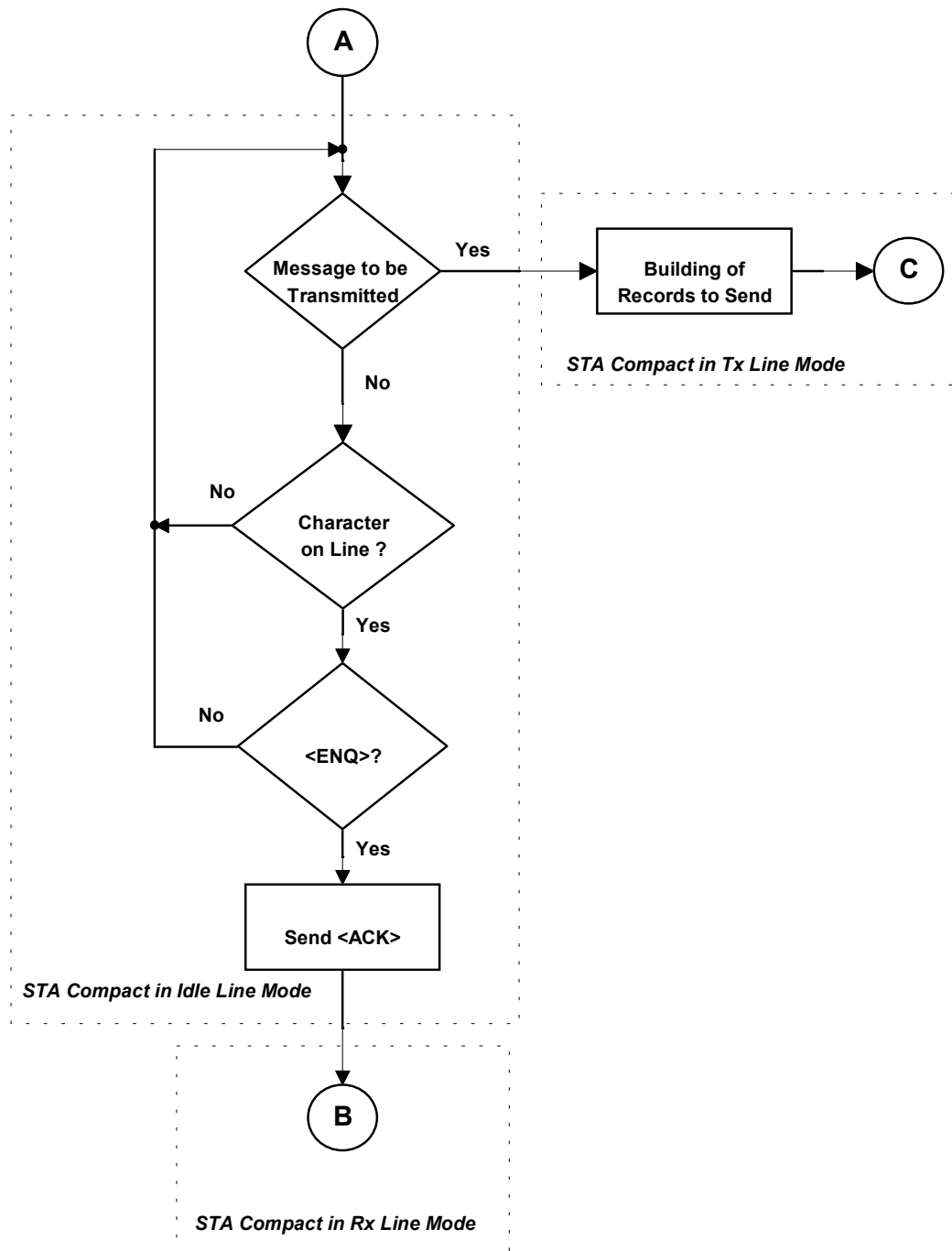


8.2. Acquisition of a Record from the Host Computer

Fig. 3 - Acquisition of a Record from the Host Computer

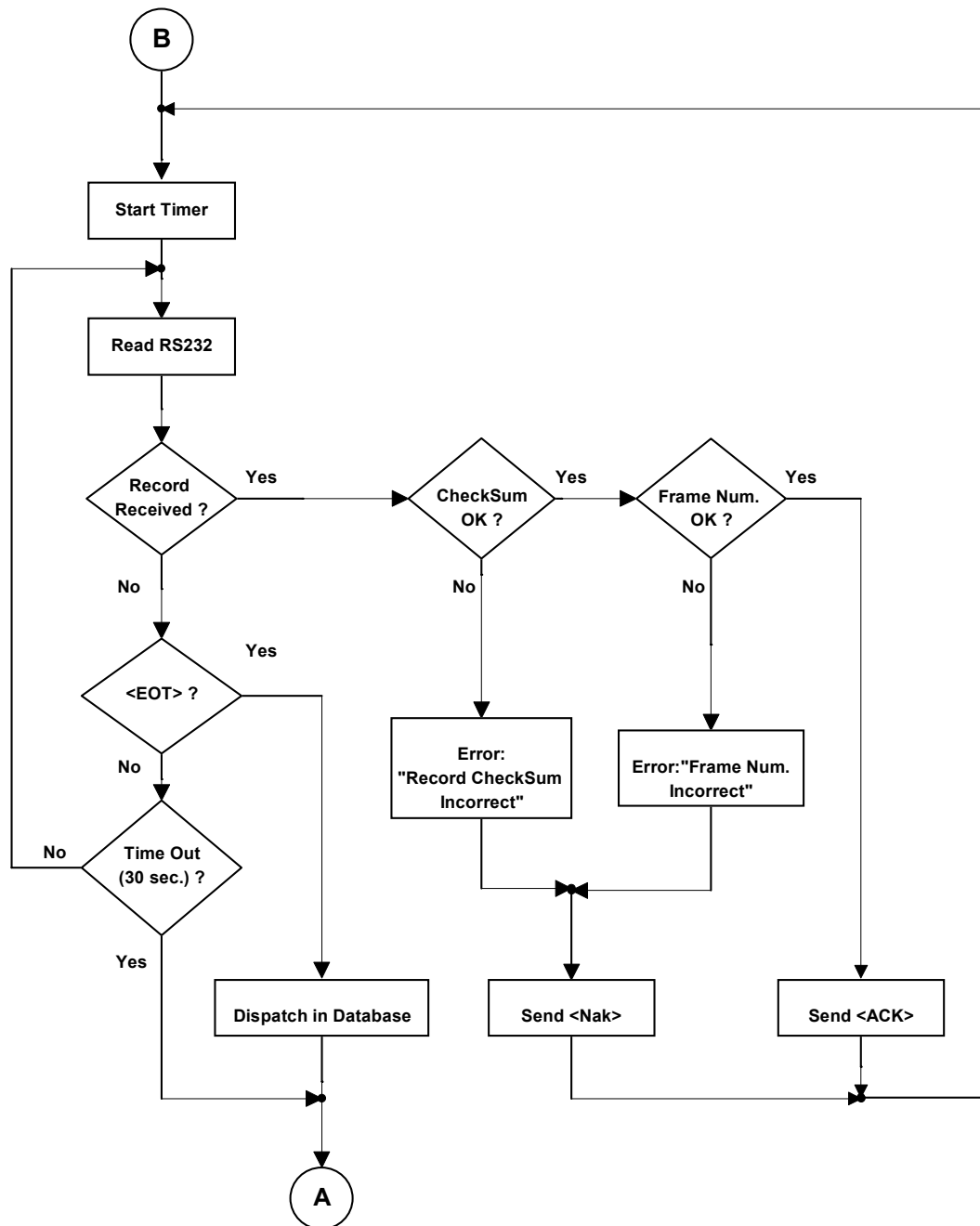

8.3. Establishment Phase (STA Compact® → Host Computer)


Fig. 4 - Establishment Phase (STA Compact® → Host Computer)

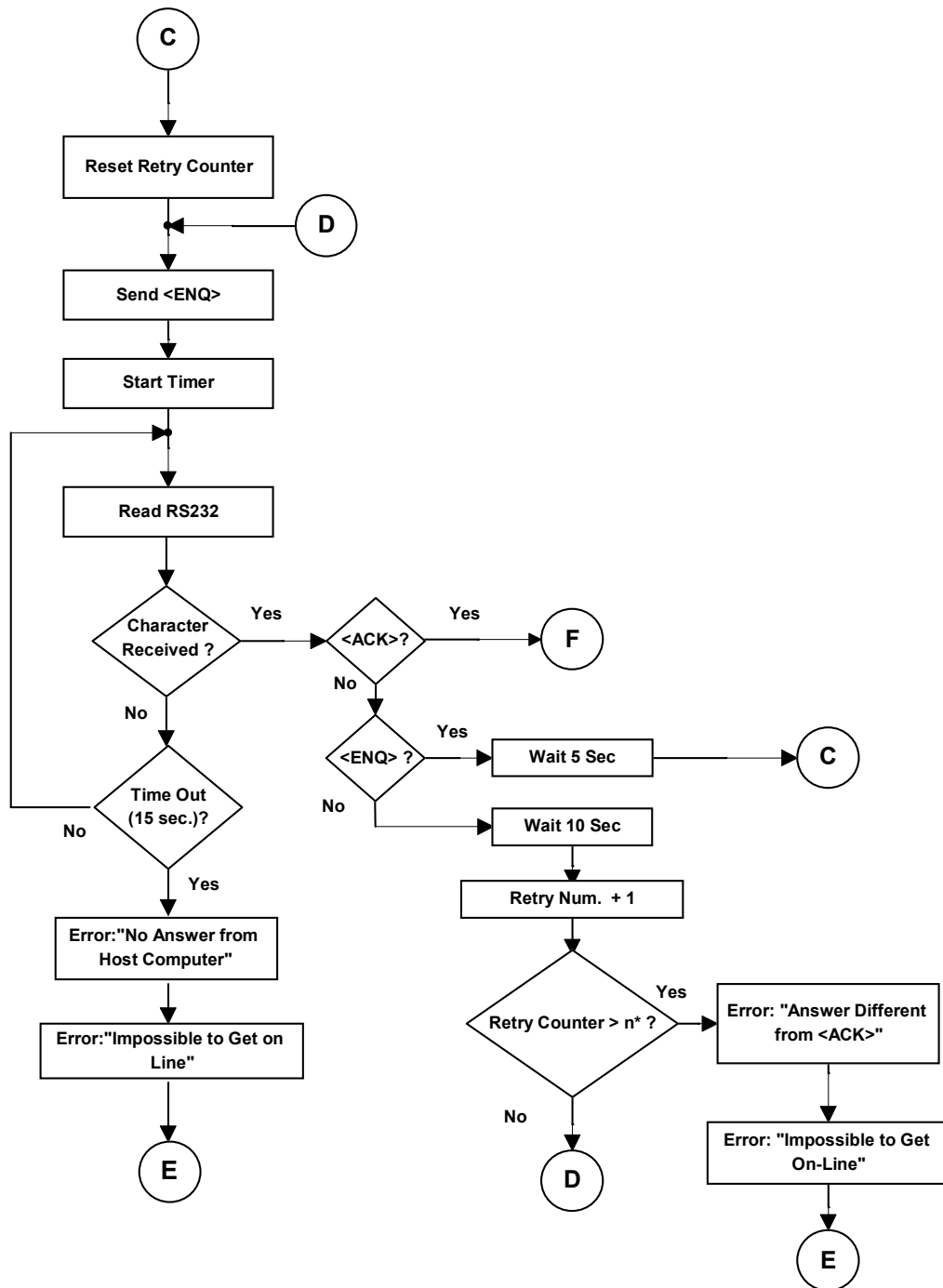

* **n** : Number of on Error Retries, defined with the **Global Options** menu, see Chapter 4. of the STA Compact® User Interface RS232-C.

8.4 Transfer Phase (STA Compact[®] → Host Computer)


Fig. 5 - Transfer Phase (STA Compact[®] → Host Computer)

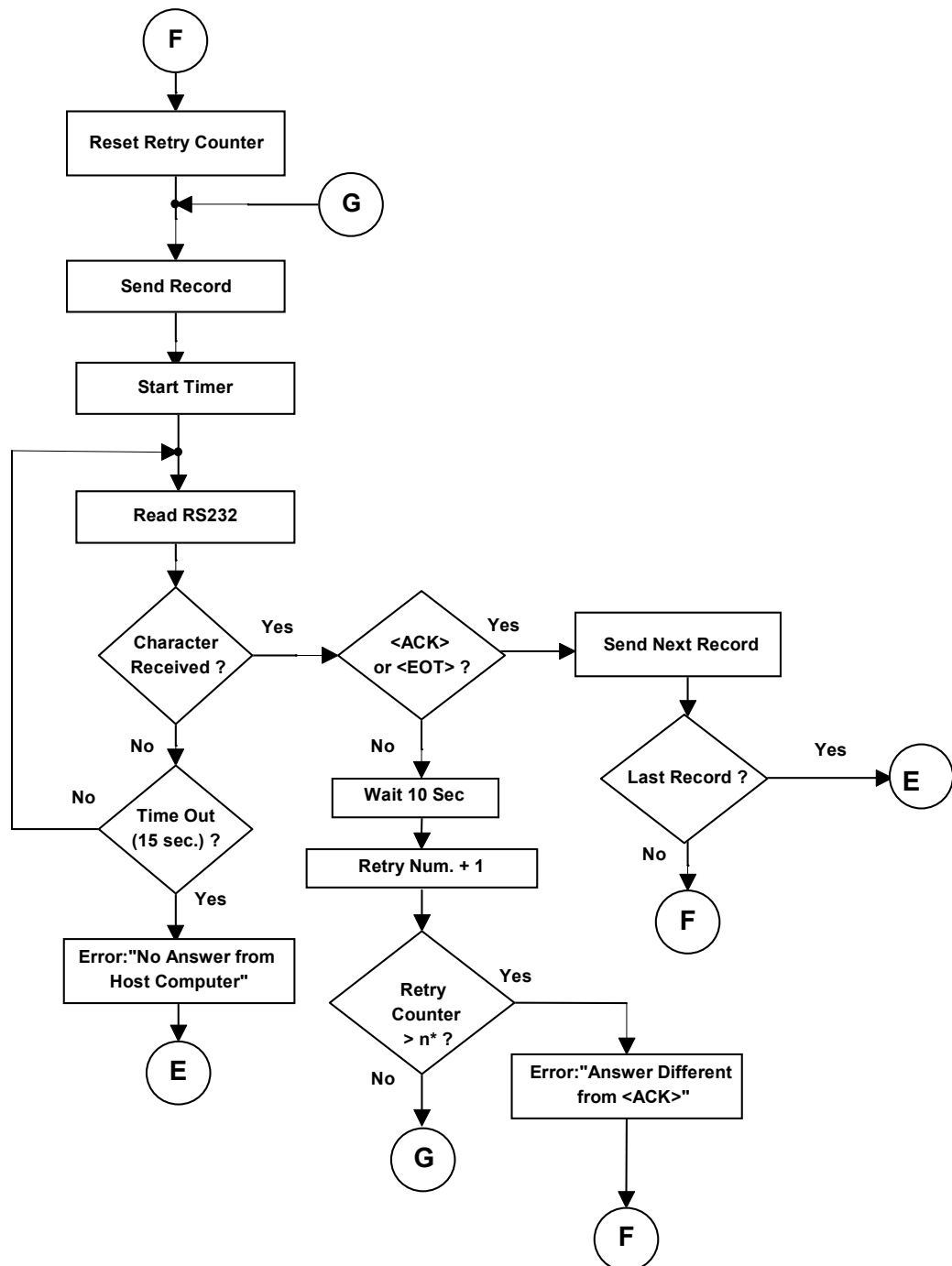

* **n** : Number of on Error Retries, defined with the **Global Options** menu, see Chapter 4.
of the STA Compact[®] User Interface RS232-C.

8.5. Termination Phase (STA Compact[®] → Host Computer)

Fig. 6 - Termination Phase (STA Compact[®] → Host Computer)

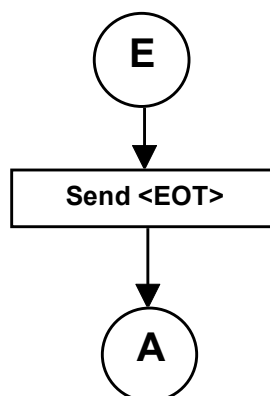

9. Examples of Traces (ASTM)

9.1. Transmission of Patient File (STA Compact[®] → Host Computer)

STA Compact [®] <ENQ>	16:09:50	STA Compact [®] requests connection.
HOST <ACK>	16:09:50	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>1H ^\& 99^2.00 P 1.00 19950227160750<CR><ETX>2A<CR><LF>	16:09:50	The STA Compact [®] sends the Header record.
HOST <ACK>	16:09:50	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>2P 1 GISCARD^Gaston^Serv.1^Gr.A <CR><ETX>5D<CR><LF>	16:09:50	STA Compact [®] sends the Patient record.
HOST <ACK>	16:09:50	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>3O 1 6 R<CR><ETX>B7<CR><LF>	16:09:50	STA Compact [®] sends the Test Order record.
HOST <ACK>	16:09:50	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>4R 1 ^1 100 % F <CR><ETX> DE<CR><LF>	16:09:50	STA Compact [®] sends the first Result record: PT = 100%.
HOST <ACK>	16:09:51	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>5M 1 A C<CR><ETX>BB<CR><LF>	16:09:51	STA Compact [®] sends the Manufacturer Information record for the first result. Error code=Validated, alarm code=Quality Control: out or range or not done.
HOST <ACK>	16:09:51	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>6R 2 ^10 10.8 sec F <CR> <ETX>5D<CR><LF>	16:09:51	STA Compact [®] sends the second Result record: PT = 10.8 seconds.
HOST <ACK>	16:09:51	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>7M 2 A C<CR><ETX>BE<CR><LF>	16:09:51	STA Compact [®] sends the Manufacturer Information record for the second result. Error code=Validated, alarm code=Quality Control: out or range or not done.
HOST <ACK>	16:09:51	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>0R 3 ^11 1.00 INR F <CR> <ETX>FF<CR><LF>	16:01:51	STA Compact [®] sends the third Result record: PT, INR=1.00.
HOST <ACK>	16:01:51	Host answers O.K. (Acknowledge).

STA Compact® <STX>1M 3 A C<CR><ETX>B9<CR><LF>	16:09:51	STA Compact® sends the Manufacturer Information record for the third result. Error code=Validated, alarm code=Quality Control: out of range or not done.
HOST <ACK>	16:09:51	Host answers O.K. (Acknowledge).
STA Compact® <STX>2R 4 ^12 12.3 Tém. F <CR> <ETX>90<CR><LF>	16:09:51	STA Compact® sends the fourth Result record: PT, Reference Time = 12.3 seconds.
HOST <ACK>	16:09:51	Host answers O.K. (Acknowledge).
STA Compact® <STX>3M 4 A C<CR><ETX>BC<CR><LF>	16:09:52	STA Compact® sends the Manufacturer Information record for the fourth result: Error code=Validated, alarm code=Quality Control: out of range or not done.
HOST <ACK>	16:09:52	Host answers O.K. (Acknowledge).
STA Compact® <STX>4R 5 ^3 4.56 g/l F <CR><ETX> FD<CR><LF>	16:09:52	STA Compact® sends the fifth Result record: Fibrinogen = 4.56 g/l.
HOST <ACK>	16:09:52	Host answers O.K. (Acknowledge).
STA Compact® <STX>5M 5 A C<CR><ETX>BF<CR><LF>	16:09:52	STA Compact® sends the Manufacturer Information record for the fifth result. Error code=Validated, alarm code= Quality Control: out of range or not done.
HOST <ACK>	16:09:52	Host answers O.K. (Acknowledge).
STA Compact® <STX>6R 6 ^30 11.9 sec F <CR> <ETX>65<CR><LF>	16:09:52	STA Compact® sends the sixth Result record: Fibrinogen = 11.9 seconds.
Site central <ACK>	16:09:52	Host answers O.K. (Acknowledge).
STA Compact® <STX>7M 6 A C<CR><ETX>C2<CR><LF>	16:09:52	STA Compact® sends the Manufacturer Information record for the sixth result. Error code=Validated, alarm code= Quality Control: out of range or not done.
HOST <ACK>	16:09:52	Host answers O.K. (Acknowledge).
STA Compact® <STX>0L 1 N<CR><ETX>03<CR><LF>	16:09:52	STA Compact® sends the Message Terminator record.
HOST <ACK>	16:09:52	Host answers O.K. (Acknowledge).
STA Compact® <EOT>	16:09:52	STA Compact® sends an END of Transmission.

9.2. Request for a Working List

STA Compact [®] <ENQ>	16:11:52	STA Compact [®] requests connection.
HOST <ACK>	16:11:52	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>1H \\& 99^2.00 P 1.00 19950227160953<CR><ETX>2F<CR><LF>	16:11:52	STA Compact [®] sends the Header record.
HOST <ACK>	16:11:52	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>2Q 1 ^ESSAI<CR><ETX>8F<CR> <LF>	16:11:52	STA Compact [®] sends the Request Information record for patient ESSAI.
HOST <ACK>	16:11:52	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>3L 1 N<CR><ETX>06<CR><LF>	16:11:53	STA Compact [®] sends the Message Terminator record.
HOST <ACK>	16:11:53	Host answers O.K. (Acknowledge).
STA Compact [®] <EOT>	16:11:53	STA Compact [®] sends an END of Transmission.

Caution : The STA Compact[®] can transmit several messages : "Request for a Working List" one after the other without allowing the intermediate return of those (see chapter 5.4. of this STA Compact[®] User Interface RS232-C and Standard E1381-91).

9.3 Working List Return

HOST <ENQ>	16:11:53	Host requests connection.
STA Compact® <ACK>	16:11:53	STA Compact® answers O.K. (Acknowledge).
HOST <STX>1H \\& 99^2.00 P 1.00 19950227161153<CR><ETX>28<CR><LF>	16:11:53	Host sends the Header record.
STA Compact® <ACK>	16:11:53	STA Compact® answers O.K. (Acknowledge).
HOST <STX>2P 1 BRUN^Didier^Essai^Site <CR><ETX>DF<CR><LF>	16:11:53	Host sends the Patient record for patient ESSAI.
STA Compact® <ACK>	16:11:54	STA Compact® answers O.K. (Acknowledge).
HOST <STX>3O 1 ESSAI ^1^2^3 R<CR> <ETX>92<CR><LF>	16:11:54	Host sends the Test Order record for patient ESSAI.
STA Compact® <ACK>	16:11:54	STA Compact® answers O.K. (Acknowledge).
HOST <STX>4L 1 N<CR><ETX>07<CR><LF>	16:11:55	Host sends the Message Terminator record.
STA Compact® <ACK>	16:11:55	STA Compact® answers O.K. (Acknowledge).
HOST <EOT>	16:11:55	Host sends an END of Transmission.

Note : This message "Working List Return" can also be used to add tests, to rerun tests or to modify the priority code.

To add tests, Test Order record has to include the local codes of tests to be added.

To rerun tests, it has to include the local codes of the tests to be rerun.

And, to modify only the priority code for a patient file already present in the STA Compact® Working File, short Test Order records can be used:

- example to modify priority code from "Routine" to "Stat" : O|1|ESSAI|||S<CR>.

9.4. Transmission of a Quality Control Result (STA Compact[®] → Host Computer)

STA Compact [®] <ENQ>	16:10:48	STA Compact [®] requests connection.
HOST <ACK>	16:10:48	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>1H ∧& 99^2.00 Q 1.00 19950227160848<CR><ETX>33<CR><LF>	16:10:48	STA Compact [®] sends the Header record with Processing ID=Q for quality control.
HOST <ACK>	16:10:48	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>2P 1 <CR><ETX>B3<CR><LF>	16:10:48	STA Compact [®] sends the Patient record.
HOST <ACK>	16:10:48	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>3O 1 12352 R<CR><ETX>7E<CR><LF>	16:10:48	STA Compact [®] sends the Test Order record for the STA [®] -COAG CONT P (12352).
HOST <ACK>	16:10:48	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>4R 1 ^1 30 % F 19950224085 100<CR><ETX>7E<CR><LF>	16:10:48	STA Compact [®] sends the Result record: PT = 30% dated 02/24/95 at 08:51.
HOST <ACK>	16:10:48	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>5M 1 A @<CR><ETX>B8<CR><LF>	16:10:48	STA Compact [®] sends the Manufacturer Information record. Error code=Validated, alarm code=none.
HOST <ACK>	16:10:48	Host answers O.K. (Acknowledge).
STA Compact [®] <STX>6L 1 N<CR><ETX>09<CR><LF>	16:10:48	STA Compact [®] sends the Message Terminator record.
HOST <ACK>	16:10:48	Host answers O.K. (Acknowledge).
STA Compact [®] <EOT>	16:10:48	STA Compact [®] sends an END of transmission.

9.5. On-Line Test

STA Compact [®] <ENQ>	16:08:14	STA Compact [®] requests connection.
HOST <ACK>	16:08:14	Host answers O.K. (Acknowledge).
STA Compact [®] <EOT>	16:08:14	STA Compact [®] sends an END of transmission.