



Procedure of Biolyte 2000 Transfer Data to Computer

Requirement:

1. RS232 Cable (Female 9-PIN and male 25-PIN),
or 1 male 25-PIN cable with an adapter.

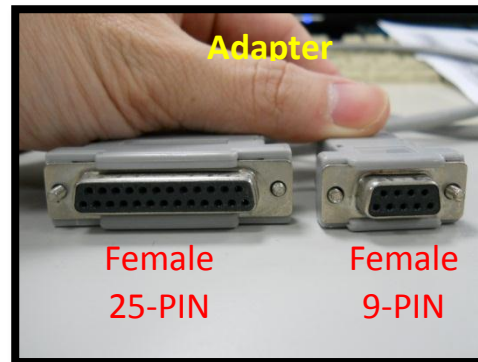

2. HyperTerminal Program.

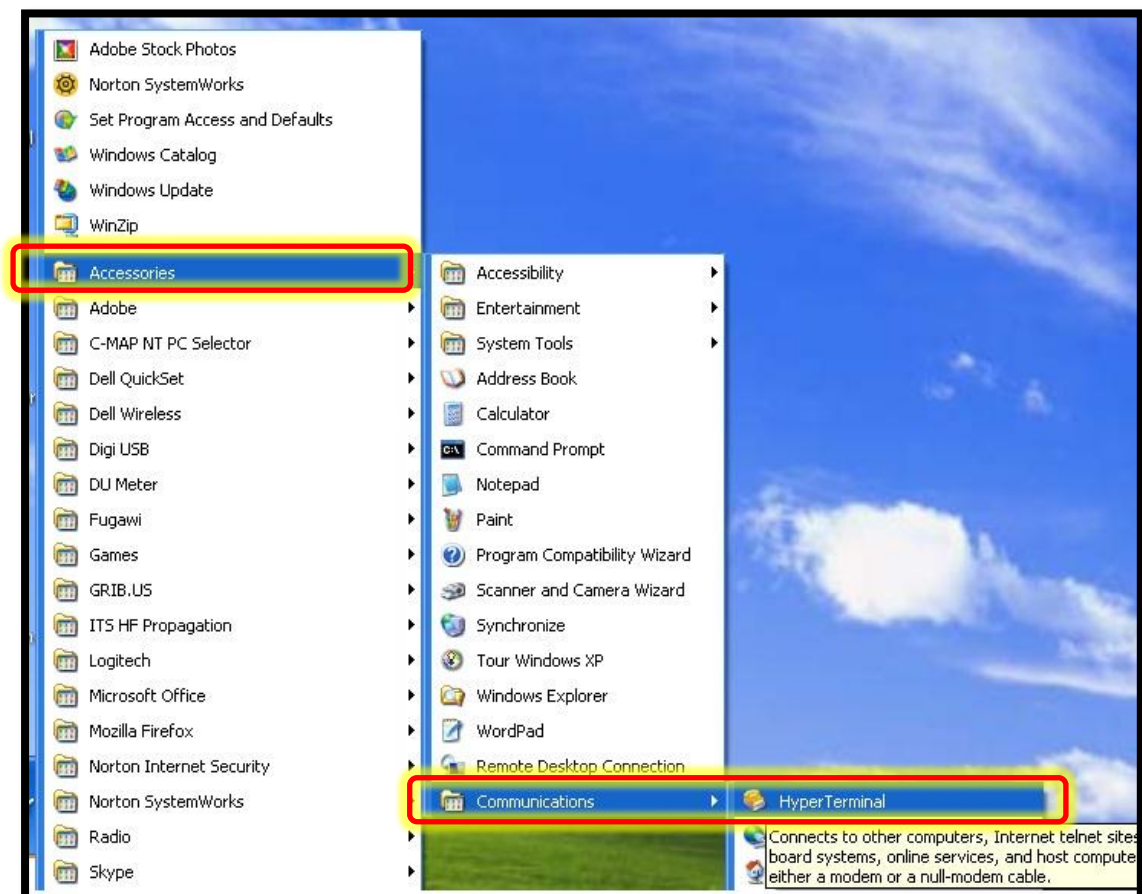


Procedure:


1. Connect Computer with Biolyte 2000 (RS232 Cable)




2. Run HyperTerminal program.

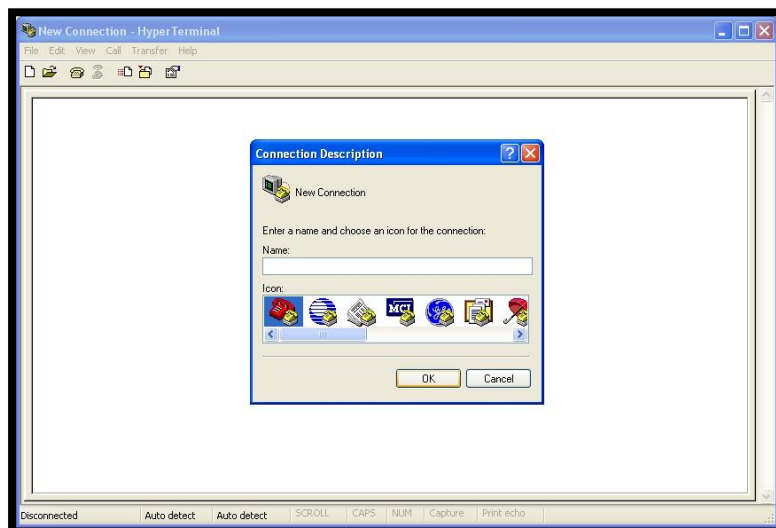

3. Choose the correct path (COM1 or COM2)

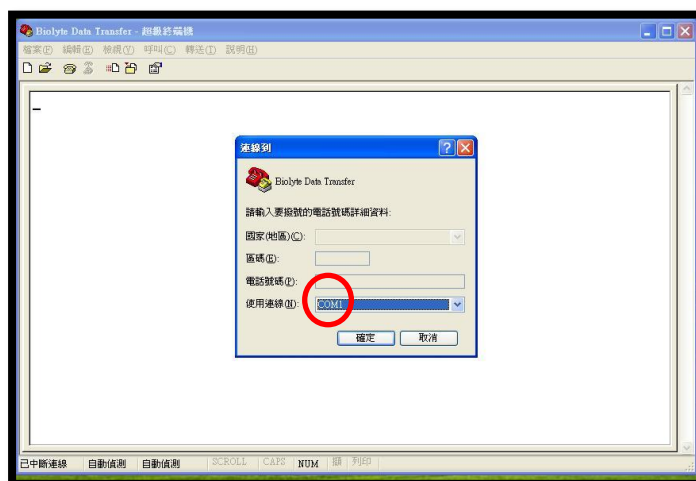

- Speed is 9600 bps.

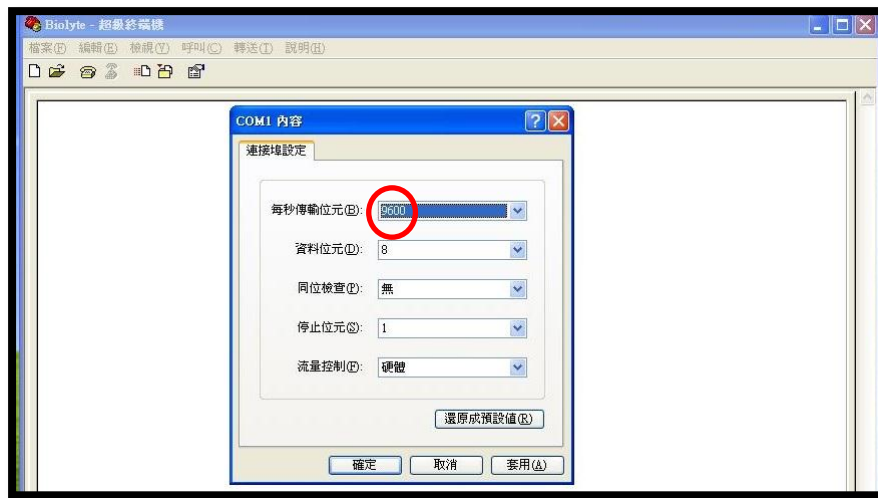

- Run Data Transfer.
(Path: Menu → 3.Data → 2.Transfer)



Then you will see this screen.



- When data are transferring press **ALT + 6 (Numerical Keyboard)** to receive each data.



7. You can find the result on your PC monitor.

